

### **500HB Installation Instructions**





PELSTAR, LLC 9500 West 55<sup>th</sup> St. McCook, IL 60525-7110 USA www.homscales.com © Pelstar, LLC 2021

weigheasier®



Thank you for purchasing this Health o meter® Professional product. To prevent patient/caregiver injury or damage to your scale, please follow the instructions in this user manual very carefully.

- The handlebars are intended for patient use only. Do not use the handlebars to transport the scale or use the scale as a handcart to transport items.
- To prevent patient injury, the patient must be attended throughout the entire weighing event. Caregivers should ensure that the patient is stable and provide support as needed when getting onto and off of the scale.
- Do not allow the patient to use the handlebars to pull themselves onto the scale.
- Excessive movement of the handlebars may affect the weight lock function.
- Do not allow patients to use excessive force to push or pull on the handlebars during the weighing process.

In no event whatsoever shall Pelstar\*, LLC be liable for damages or injuries whatsoever arising out of or connected with the assembly, use, or misuse of its products.

**Tools Included:** ¼" Drill Bit\*, Allen Wrench, Ratchet Allen Wrench

### **Tools Required:** Utility Knife, Power Drill\*

(\*Note: Power drill and drill bit are not required when mounting handlebars onto scales with factory installed mounting holes and cover caps. See photo on page 4.)

### **Parts List:**













1/4-20 X 2 1/4 Button Head Screws, 8 Each



1/4-20 X 1/2 Button Head Screws, 4 Each



500HB Mount Bracket, 2 Each








	Gently place the fully assembled scale on a flat, level surface on its side.	
2.	Remove the 2 Rear Screws (closest to the wheels) that secure the Scale Column to the Base.	
3.	Loosen but <b>DO NOT REMOVE</b> the 2 Forward Screws (closest to the scale base). Loosen Screws approximately ¾". Note, the Load Cell Cable does not need to be removed.	Loosen Screws
4.	Loosen Screws approximately ¾".	Loosen Screws ¾"
5.	Tilt the front edge of the Base toward the top of the Scale opening a gap at the rear of the Scale.  Caution: Ensure the load cell cable does not disconnect from the load cell port.	Gap at Rear of Base  Gap at Rear of Base
6.	Through the opening in the gap, slide the Spacer into the Gap. Position the Spacer so that the open end of the Spacer slides around the Load Cell Cable.	

7. Slide Spacer into Base and position Spacer so that Mount Holes in Spacer align with Holes in Base



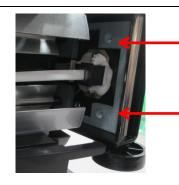


8. Insert the two Screws removed in step 3 into the Mount Holes on bottom of Base. Then tighten all 4 Screws securely.






Note: If your scale has the factory installed mounting holes and cover caps as shown, follow installation step 9A and proceed to step 17.

# For scales without cover caps complete step 9.B then proceed to step 10.

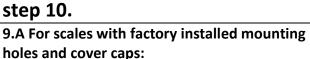

Place the drill guide onto the bottom of the scale platform as shown. Obtain the included Allen wrench and insert the wrench into the holes in the drill guide to push the cover caps upward. Remove all (4) caps from this side of the scale and discard. Repeat these steps for the other side of the scale. Remove and discard the drill guide then skip to step 17.

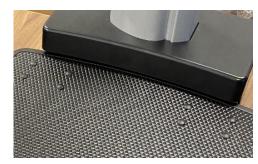


Scales with cover caps, perform step 9.A then proceed to step 17.

Align edge of drill guide with edge of scale platform

## 9.B For scales <u>without</u> factory installed mounting holes and cover caps:

Place the drill guide onto the bottom of the scale platform as shown. Position drill guide so that the arrow on the guide points to center of the scale. Align edge of drill guide with edge of platform.

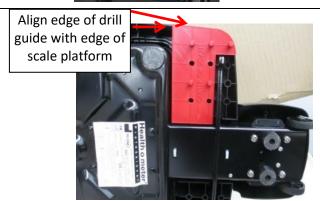


10. There are posts on both sides of the drill guide. These posts must slide into holes in platform.	
11. Drill guide must be positioned so that posts slide <b>fully</b> into holes on bottom of scale platform.	
12. With drill guide properly positioned on bottom of scale platform, use a power drill and the provided ¼" drill bit to drill 4 holes through the bottom of the scale platform.	
13. Flip the drill guide (drill guide is used on both sides of the scale platform) and position it on the other side of the bottom of the platform.	To mage
14. With drill guide properly positioned on bottom of scale platform, use a power drill and the provided ¼" drill bit to drill 4 holes through the other of the side bottom of the scale platform.	

[	
15. Remove and discard drill guide. Carefully place the scale upright on the floor. Platform will now have 8 holes for mounting the handlebars.	
16. Use a utility knife or similar tool to remove plastic burrs from the top of the scale platform.	
17. Using the provided alcohol wipe, clean the entire top surface of the scale platform.	
18. Obtain the handlebar and base. Align the bottom of the handlebar with the channels on the base. Slide the handlebar onto the channel until the bottom of the handlebar rests on the base. Ensure the holes in the handlebar align with the holes in the channel.	
19. Obtain the four ¼-20 x ½" Button Head Screws. 20. Insert all four screws into the back of the handlebar and into the channel until threads are engaged. Use the included ratchet allen wrench to fully tighten and secure the screws. Note:  Only tighten once all four screws have been inserted into the threads. Ensure the switch on the ratchet allen wrench is set to the tightening action.	
20. Place the handlebar and base assembly around the display head as shown and lower onto the base of the scale. Position the base onto the scale platform by aligning the holes in the plate with the holes in the drilled in the platform.  Note: There is a strip of double-sided tape on the bottom of the handlebar base. This will be removed in step 22.	

21. Insert two of the provided 1/4-20 X 2 1/4" screws into the plate on the handlebars and into the scale platform, one screw on the outside of both sides of the handlebar channel. Insert the screws manually or using the included ratchet allen wrench. Do not fully tighten the screw at this stage.	
22. Locate the strip of double-sided tape on the bottom of the plate.	
23. Slightly tilt the base up to allow for the removal of the tape backing strip. Peel off the backing strip under the plate. Lower the plate/handlebars onto the scale platform and apply light pressure for 30 seconds.	
24. Using the included ratchet allen wrench, insert the remaining six button head screws through the base and into the scale platform. Do not completely tighten.	
25. Access the bottom of the platform by placing the scale carefully on its side or up on a table.  Obtain one of the silver mounting brackets and position it so that the flat side of the bracket faces the bottom of the platform.	
26. Position the bracket as shown and thread the screw into the threaded nut on the bracket.	

27. **Note:** When mounting the brackets underneath the platform, ensure that the silver brackets are not touching the black stabilizer bar. If necessary, loosen the screws to adjust, then retighten.

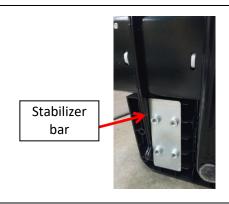

28. Repeat the process with the second mounting bracket. Then tighten all 8 screws completely. Assembly is complete.



### **Limited Warranty**

### What does the Warranty Cover?

This Health o meter® Professional product is warranted from date of purchase against defects of materials or in workmanship for a period of two (2) years. If product fails to function properly, return the product, freight prepaid and properly packed to Pelstar, LLC (see "To Get Warranty Service", below, for instructions). If the manufacturer determines that a defect of material or in workmanship exists, the customer's sole remedy will be replacement of the product at no charge. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. All replaced parts are covered only for the original warranty period.

### Who is Covered?


The original purchaser of the product must have proof of purchase to receive warranty service. Please save your invoice or receipt. Pelstar dealers or retail stores selling Pelstar products do not have the right to alter, or modify or any way change the terms and conditions of this warranty.

#### What is Excluded?

Your warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, abuse including tampering, damage in transit, or unauthorized repair or alternations. Further, the warranty does not cover natural disasters, such as fire, flood, hurricanes and tornadoes. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country, state to state, province to province or jurisdiction to jurisdiction.

To get Warranty Service make sure you keep your sales receipt or document showing proof of purchase. Call (+1) 800-638-3722 or (+1) 708-377-0600 to receive a return authorization (RA) number, which must be included on the return label. Attach your proof of purchase to your defective product along with your name, address, daytime telephone number and description of the problem. Carefully package the product and send with shipping and insurance prepaid to:

Pelstar, LLC
Attention R/A#\_\_\_\_
Return Department
9500 West 55<sup>th</sup> Street
McCook, IL 60525

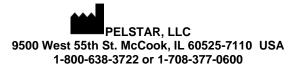

PLEASE REGISTER YOUR PRODUCT FOR WARRANTY COVERAGE AT:

www.homscales.com

Health o meter<sup>®</sup> is a registered trademark of Sunbeam Products, Inc. used under license. Health o meter<sup>®</sup> Professional products are manufactured, designed, and owned by Pelstar, LLC. We reserve the right to improve, enhance, or modify Health o meter<sup>®</sup> Professional product features or specifications without notice.