

## THE 2600KL


# DIGITAL WHEELCHAIR SINGLE RAMP SCALE

## THE 2610KL

# DIGITAL WHEELCHAIR DUAL RAMP SCALE

The all-digital 2600KL is our most popular wheelchair scale that supports up to 1000 pounds and provides the optimum in quality, reliability, and ease of use. The platform folds for convenient storage and the casters (the only scale in the market that has them!) provide easy mobility. The wide and low profile platform accommodates wheelchairs and provides a comfortable weighing experience for patients. For easy on/off access the 2610KL offers dual ramps. Designed to meet the needs of doctor's offices, hospitals, clinics, and acute and long-term care facilities, the scale's keypad features antimicrobial protection, a display head that swivels 180°, a large non-skid low-profile platform, and wheels that enable portability.

The scale's high-contrast color screen provides a user friendly interface that includes BMI, hold/release, reweigh, recall, zero, KG/LB, tare, pre-tare and more. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently locked on the scale. The scale provides easy integration and seamless communication with today's EMR applications, the 2600KL and 2610KL include a USB interface for connectivity with a computer, EMR software, or other electronic devices via various Health o meter® Professional Connectivity Solutions.

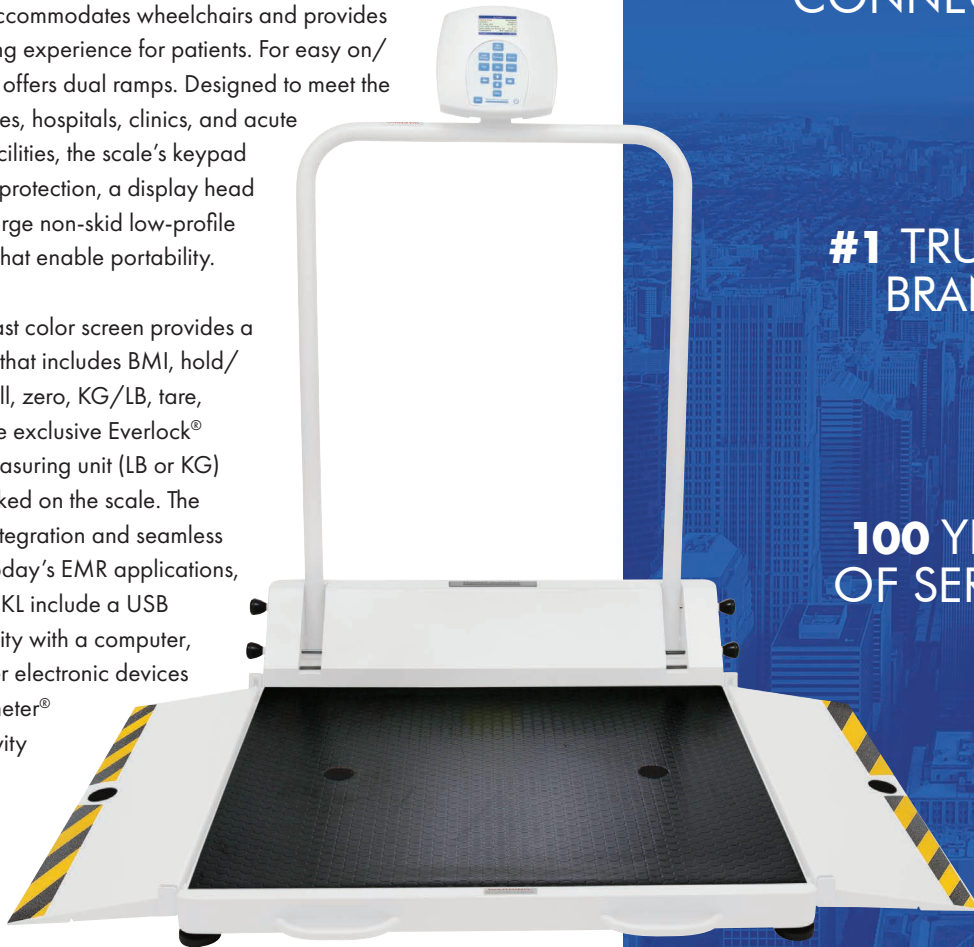

2610KL

## THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 2600KL is our most popular wheelchair scale.

**120+** LOCAL  
FIELD SPECIALISTS

**100%**  
CONNECTIVITY

**#1** TRUSTED  
BRAND

**100** YEARS  
OF SERVICE

**Health o meter®**  
Professional

Health o meter® Professional 2600KL & 2610KL

weigh**easier**®

# TECHNICAL SPECIFICATIONS

## Capacity

1000 lb / 454 kg

## Resolution

0.2 lb / 0.1 kg

## EMR connectivity

Via USB, optional Pelstar® wireless technology

## Display

Digit Size: 3/8" / 10 mm  
Screen Size: 2 3/4" x 1 1/2" High-Contrast  
Color TFT-LCD Screen

## Power Source

100-240V adapter (ADPT30) included, 6 D batteries not included

## Functions

LB/KG conversion  
LB/KG lock out  
BMI  
Zero  
Tare  
Pre-tare  
Hold/release  
Reweigh  
Recall  
Auto zero  
Auto off  
Time/date  
Variable auto off time  
Audible/mute sound option  
Parent/child weighing function  
Everlock®

## Platform Size (w x d x h)

36" x 32 1/4" x 2 1/2"  
(914 mm x 819 mm x 64 mm)

## Ramp Size

32 1/4" x 9 7/8"  
(819 mm x 248 mm)

## 2600KL Product Footprint

48 1/8" x 43 7/8" x 49 3/8"  
(1222 mm x 1114 mm x 1254 mm)

## 2610KL Product Footprint

56 3/4" x 43 7/8" x 49 3/8"  
(1441 mm x 1114 mm x 1254 mm)

## Product Weight

2600KL: 93 lb / 42 kg, 2610KL: 98 lb / 44 kg

## Portability

4 wheels (2 casters)  
for ease of movement

## Warranty

2-year limited manufacturer warranty

# ORDERING INFORMATION

MODEL	DESCRIPTION
2600KL	Digital wheelchair single ramp scale, LB & KG
2610KL	Digital wheelchair dual ramp scale, LB & KG*
2620KL	Digital Wheelchair Scale with Dual Ramps and Extra-Large Platform, ADA Compliant, LB & KG*
2600KG	Digital wheelchair single ramp scale, KG only
2610KG	Digital wheelchair dual ramp scale, KG only*
2620KG	Digital Wheelchair Scale with Dual Ramps and Extra-Large Platform, ADA Compliant, KG only*
2600KL-BT	2600KL with built-in Pelstar® wireless technology
2610KL-BT	2610KL with built-in Pelstar® wireless technology*
2600KG-BT	2600KG with built-in Pelstar® wireless technology
2610KG-BT	2610KG with built-in Pelstar® wireless technology*
2600ADAKIT	ADA Compliance Kit for Gray 2600 and 2610 Series Scales*
2600ADAKIT-WHITE	ADA Compliance Kit for White 2600 and 2610 Series Scales*
2622ADAKIT	ADA Wheelchair & Stand-on Compliance Kit for White 2600 and 2610 Series Scales*
B2600RAMP	Optional secondary ramp for White 2600 Series Scales
43-00063	Optional secondary ramp for Gray 2600 Series Scales
SS-2600KL	ScaleSurance® 2 year extended warranty for 2600 series*
SS-2610KL	ScaleSurance® 2 year extended warranty for 2610 series*
SS-2620KL	ScaleSurance® 2 year extended warranty for 2620 series*
ADPT30	Adapter included, item # for replacement or spare
PLUGSET30	Plug adapters for use outside North America
BTPRINTERKIT-01	Wireless thermal printer kit

Scale also includes connectivity capabilities. To learn about connectivity options, visit [hom scales.com/ConnectivitySolutions](http://hom scales.com/ConnectivitySolutions)

To learn more, call **1-800-815-6615**  
or visit [hom scales.com/2600KL](http://hom scales.com/2600KL)  
[hom scales.com/2610KL](http://hom scales.com/2610KL)



2600KL



2610KL folded

**Health o meter®**  
Professional



\*Not available in all countries.

2600KL\_2610KL\_F&B\_20260604