


THE 3001KL-AM

ANTIMICROBIAL DIGITAL PLATFORM SCALE

The all-digital 3001KL-AM supports up to 1000 pounds and provides the optimum quality, reliability, and ease of use designed to meet the needs of doctor's offices, hospitals, clinics, and acute and long-term care facilities. The scale features a large non-skid low profile platform, "live" handrails for patient stability, and rugged scuff proof wheels for easy repositioning. The scale also offers an antimicrobial keypad, handrails, and transport handle to give patients and providers additional safety.

The scale's high-contrast color screen provides a user-friendly interface that includes BMI, hold/release, reweigh, recall, zero, KG/LB, tare, pre-tare, and more. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently locked on the scale.

Antimicrobial
Keypad
Handrails
Transport Handle

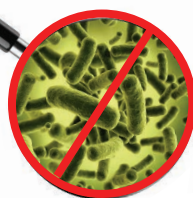

The 3001KL-AM offers two optional height rod options with quick and easy to assemble brackets, a digital height rod 245EHR-3001 or a mechanical height rod 3001HR. The 245EHR-3001 market-leading digital height rod provides the most accurate height measurement in its class, with a graduation/resolution of 1mm/0.1" and a full measuring range of 11 3/4" – 79 1/2" / 30cm – 202cm and a digital range of 3ft 11.2" – 6ft 7.5" / 120 cm - 202 cm. The 3001HR features a durable height rod construction and an easy-slide headpiece that allows caregivers to quickly and effortlessly measure and read a patient's height. The 3001HR provides a measuring range of 2"-90" in 1/16" increments (6 cm-230cm in 1 mm increments).

The scale provides easy integration and seamless communication with today's EMR applications. The 3001KL-AM includes a USB interface for connectivity with a computer, EMR software, or other electronic devices via various Health o meter® Professional Connectivity Solutions.

THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 3001KL-AM is your choice for safety and reliability.

120+ LOCAL
FIELD SPECIALISTS

100%
CONNECTIVITY

#1 TRUSTED
BRAND

100 YEARS
OF SERVICE



Health o meter®
Professional

Health o meter® Professional 3001KL-AM

weigh**easier**®

TECHNICAL SPECIFICATIONS

Capacity

1000 lb / 454 kg

Resolution

0.2 lb / 0.1 kg

EMR Connectivity

Via USB, optional Pelstar® wireless technology

Display

Digit Size: 3/8" / 10 mm

Screen Size: 2 3/4" x 1 1/2" High-Contrast

Color TFT-LCD Screen

Power Source

100-240V adapter (ADPT30)

included, 6 D batteries not included

Functions

LB/KG conversion

LB/KG lock out

BMI

Zero

Tare

Pre-tare

Hold/release

Reweigh

Recall

Auto zero

Auto off

Time/date

Variable auto off time

Audible/mute sound option

Parent/child weighing function

Everlock® - permanent weight lock

Platform Size (w x d x h)

22" x 15 3/4" x 2 5/8"

(559 mm x 400 mm x 64 mm)

Product Footprint (w x d x h)

22" x 33 1/8" x 50"

(559 mm x 841 mm x 1270 mm)

Product Weight

46 lb / 21 kg

Portability

2 wheels for ease of movement

Warranty

2-year limited manufacturer warranty

ORDERING INFORMATION

MODEL	DESCRIPTION
3001KL-AM	Antimicrobial digital platform scale, LB & KG
3001KG-AM	Antimicrobial digital platform scale, KG only
3001KL-AM-BT	3001KL-AM with built-in Pelstar® wireless technology
3001KG-AM-BT	3001KG-AM with built-in Pelstar® wireless technology
3001KL-AMHR	3001KL-AM with mechanical height rod
3001KG-AMHR	3001KG-AM with mechanical height rod
3001KL-AMHR-BT	3001KL-AMHR with built-in Pelstar wireless technology
245EHR-3001	Optional digital height rod for 3001 models*
3001HR	Optional mechanical height rod for 3001 models
SS-3001KL-AM	ScaleSurance 2 year extended warranty (3001 models)**
ADPT30	Adapter included, item # for replacement or spare
PLUGSET30	Plug adapters for use outside North America
BTPRINTERKIT-01	Wireless thermal printer kit

Scale also includes connectivity capabilities. To learn about connectivity options, visit homscales.com/ConnectivitySolutions

Optional Height Rods Available



To learn more, call **1-800-815-6615**
or visit homscales.com/3001KL-AM

Health o meter®
Professional



* Compatible with 3001 series scales including a date code with the last two digits of 21 or higher. An upgrade option is available for models including a date code with the last two digits of 17 or higher. Not available for AMX series scales. For assistance, contact Customer Service at 800-815-6615.

** Not available in all countries