Digital Pediatric Tray Scale

B O S S Φ



This easy-to-use scale is perfect for busy pediatric offices or clinics where infants are seen regularly. The tray has a convenient measuring tape to determine the baby's length during the weighing session.

Advantages

EMR Connectivity

Scale has the ability to interface with a computer, EMR software or other electronic device via various Health o meter® Professional Connectivity Solutions

Measuring Tape Included

Easy to read measuring tape allows user to measure both weight and height simultaneously

Motion Detection

Scale has the ability to weigh accurately despite movement of the subject on the scale

Display Head 17" Above Scale Base Display head is positioned 17" above the scale base to allow for an unobstructed view of the weight measurement

Value-Added Features

Operating functions (LB/KG Conversion, LB/KG Lock Out, Zero Out/Tare, Hold/ Release, Auto Zero, Auto Off) included at no additional cost



Digital Pediatric Tray Scale



Health o meter® Professional Scales 9500 West 55th Street USA McCook, IL 60525-7110 Phone: 1-800-815-6615 Fax: 1-800-654-7330




Specifications

EMR

MET

- Capacity: 50 lb / 23 kg
- Resolution: 0 20 lb 0.2 oz / 20 50 lb / 0.5 oz; 0 9 kg / 5g; 9 23 kg / 10 g
- Connectivity via USB, optional Pelstar wireless technology
- 1" / 25 mm LCD display
- Functions: LB/KG Conversion, LB/KG Lock Out, Zero, Tare, Hold/Release, Auto Zero, Auto Off
- Tray Size (w x d x h): 24 ¹/₈" x 14 ¹/₂" x 2 ⁵/₈"
 613 mm x 368 mm x 67 mm
- Product Footprint (w x d x h): 24 ¹/₈" x 21" x 23 ¹/₈"
 - 613 mm x 533 mm x 587 mm
- Product Weight: 20 lb / 9 kg
- Measuring Tape: 0" 23" / 0 cm 58 cm
- 6 AA batteries included, Optional adapter (order # ADPT31)
- 2 year limited warranty

Model	Description
522KL	Digital pediatric tray scale, pounds & kilograms
522KG	Digital pediatric tray scale, kilograms only
522KL-HR	522KL with attachable mechanical height rod
522KL-EHR	522KL with attachable digital height rod
SS-522KL	ScaleSurance 2 year extended warranty for 522KL & 522KG*
522HR	Attachable mechanical height rod for 522KL
522EHR	Attachable digital height rod for 522KL
ADPT31	Optional adapter
PELSTARONE-P	Wireless thermal printer

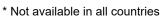



522KL_F&B_20200901