

## THE 2900KL-AM

# ANTIMICROBIAL WALL-MOUNTED FOLD-UP WHEELCHAIR SCALE

The number-one-selling medical scale company in North America offers an antimicrobial wall-mounted fold-up wheelchair scale that supports up to 1000 pounds and provides the optimum in quality, reliability, and ease of use. The low profile scale extends off the wall 3.5" when folded up and is excellent for facilities that need to save space. The large easy-to-lift platform accommodates most wheelchairs and provides a comfortable weighing experience for patients. The dual ramps offer easy on, easy off access and the anti-slip caution tape delivers enhanced safety. Other safety features include antimicrobial protection, a raised edge, a recessed handle, and no pinch points, which provide both the clinician and patient confidence during the weighing process.

Designed to meet the needs of hospitals, clinics, acute, long-term care facilities, doctor's offices, and anywhere space is limited. The scale gives providers flexibility with a wall-mounted display that can be mounted on either side of the scale. The 2900KL-AM's high-contrast color screen provides a user-friendly interface that includes BMI, hold/ release, reweigh, recall, zero, KG/LB, tare, pretare and more. The exclusive Everlock® function allows the measuring unit (LB or KG) to be permanently locked on the scale.

The scale provides easy integration and seamless communication with today's EMR applications. The 2900KL-AM includes a USB interface for connectivity with a computer, EMR software, or other electronic devices via various Health o meter® Professional Connectivity Solutions.

**Health o meter®**  
Professional

Health o meter® Professional 2900KL-AM

## THERE'S STRENGTH IN OUR NUMBERS.


When you add up the numbers, it's not surprising that the Health o meter® Professional Scales is the industry-leading name in medical scales.

**120+** LOCAL  
FIELD SPECIALISTS

**100%**  
CONNECTIVITY

**#1** TRUSTED  
BRAND

**100** YEARS  
OF SERVICE

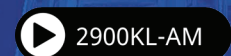

weigh**easier**®



# TECHNICAL SPECIFICATIONS

## Capacity

1000 lb / 454 kg

## Resolution

0.2 lb / 0.1 kg

## EMR connectivity

Via USB, optional Pelstar® wireless technology

## Display

Digit Size: 3/8" / 10 mm  
Screen Size: 2 3/4" x 1 1/2" High-Contrast  
Color TFT-LCD Screen

## Power Source

100-240V adapter (ADPT30) included, 6 D batteries not included

## Functions

LB/KG conversion  
LB/KG lock out  
BMI  
Zero  
Tare  
Pre-tare  
Hold/release  
Reweigh  
Recall  
Auto zero  
Auto off  
Time/date  
Variable auto off time  
Audible/mute sound option  
Parent/child weighing function  
Everlock® - permanent unit lock

## Platform Size (w x d x h)

29 1/2" x 32 1/2" x 2 3/4"  
(749 mm x 819 mm x 70 mm)

## Ramp Size (w x d)

6 1/4" x 34"  
(159 mm x 864 mm)

## Product Footprint Unfolded (w x d x h)

58" x 37" x 37"  
(1473 mm x 940 mm x 940)

## Product Footprint Folded (w x d x h)

58" x 3 1/2" x 37"  
(1473 mm x 89 mm x 940 mm)

## Product Weight

91 lb / 41 kg

## Warranty

2-year limited manufacturer warranty

# ORDERING INFORMATION

MODEL	DESCRIPTION
2900KL-AM	Antimicrobial wall-mounted fold-up wheelchair scale, LB & KG
2900KG-AM	Antimicrobial wall-mounted fold-up wheelchair scale, KG only
2900KL-AM-BT	2900KL-AM with built-in Pelstar® wireless technology
2900KG-AM-BT	2900KG-AM with built-in Pelstar® wireless technology
2920KL-AM	Antimicrobial wall-mounted fold-up wheelchair scale, ADA Compliant, LB & KG
2920KG-AM	Antimicrobial wall-mounted fold-up wheelchair scale, ADA Compliant, KG only
SS-2900KL-AM	ScaleSurance® 2 year extended warranty for 2900 series*
SS-2920KL-AM	ScaleSurance® 2 year extended warranty for 2920 series*
ADPT30	Adapter included, item # for replacement or spare
PLUGSET30	Plug adapters for use outside North America
BTPRINTERKIT-01	Wireless thermal printer kit

Scale also includes connectivity capabilities. To learn about connectivity options, visit [homscapes.com/ConnectivitySolutions](http://homscapes.com/ConnectivitySolutions)



**Easy to Lift Platform with a Recessed Handle**

**Health o meter®**  
Professional



\*Not available in all countries.

2900KL-AM\_F&B\_20251201