


THE 2210KL-AM-C

ANTIMICROBIAL HIGH-RESOLUTION DIGITAL NEONATAL/PEDIATRIC TRAY SCALE

The 2210KL-AM-C neonatal scale was designed for use in the NICU, nursery, delivery room, and pediatric practices, where a high degree of accuracy is critical to the smallest of patients. The scale gives providers and patients safety with an Antimicrobial body, tray, handles, and keypad. The smooth-edged removable tray has a patent-pending mounting system that eliminates breakage, making it easy to clean. The built-in measuring tape allows users to measure both weight and height simultaneously. The scale will quickly lock despite the movement of the infant. The backlit display makes reading a patient's weight easier at almost any angle and light condition.

The high 1 gram resolution scale features a capacity of 20,000 grams/20 kg/45lb. The 2210KL-AM-C's user-friendly digital interface includes breast milk intake, auto-weight hold, hold/release, zero, tare, pre-tare, reweigh, recall, GM/KG/LB conversion and lockout, auto-zero, and auto-off. The

scale's motion detection gives it the ability to weigh accurately despite the movement of the patient. The exclusive Everlock® function allows the measuring units (LB, KG, or G) to be permanently locked on the scale.

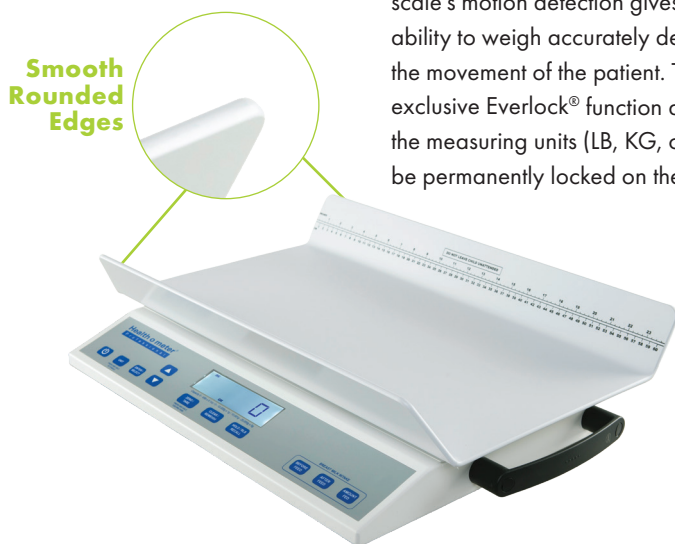

The scale provides easy integration and seamless communication with today's EMR applications. The 2210KL-AM-C includes a USB interface for connectivity with a computer, EMR software, or other electronic devices via various Health o meter® Professional Connectivity Solutions. Options such as kilograms only, grams only, and built-in wireless models are available.

Health o meter®
Professional

Health o meter® Professional 2210KL-AM-C

THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 2210KL-AM-C is your choice for accuracy and reliability

**1 GRAM
RESOLUTION**

**120+ LOCAL
FIELD SPECIALISTS**

**100%
CONNECTIVITY**

**#1 TRUSTED
BRAND**

**100 YEARS
OF SERVICE**

weigh**easier**®

TECHNICAL SPECIFICATIONS

Capacity

45 lb / 20 kg / 20,000 gm

Resolution

1 g (to 10,000 g);

5 g (10,001 to 20,000 g)

EMR Connectivity

Via USB, optional Pelstar® wireless technology

Backlit LCD Display

5 1/4" x 1 " (133 mm x 35 mm)

4 levels of brightness

Batteries

6 D batteries included

Functions

LB / KG / GM conversion

LB / KG / GM lock out

Zero

Tare

Pre-tare

Reweigh

Recall

Auto zero

Auto off

Motion detection

Everlock®

Hold / release

Variable auto-off Time

Audible/mute sound option

Disable/enable backlight for display

Enabling unit view with unit of measure Locked

Auto-weight hold for 30 to 60 seconds

Breast milk intake

Tray Size (w x d x h)

24" x 15 1/4" x 2 3/4"

(610 mm x 388 mm x 71 mm)

Product Footprint

24 3/4" x 17 3/4" x 6"

(630 mm x 450 mm x 153 mm)

Product Weight

28 lb / 13 kg

Measuring Tape

0" - 24" / 0 cm - 60 cm

Portability

Heavy duty handles for ease of moving

2210CART - optional rolling cart

Warranty

2-year limited manufacturer warranty

ORDERING INFORMATION

| MODEL | DESCRIPTION |
|----------------|--|
| 2210KL-AM-C | Antimicrobial high resolution digital neonatal / pediatric tray scale, LB & KG |
| 2210KG-AM-C | Antimicrobial high resolution digital neonatal / pediatric tray scale, KG only |
| 2210KL-AM-BT-C | 2210KL-AM-C with built-in Pelstar® wireless technology |
| 2210KG-AM-BT-C | 2210KG-AM-C with built-in Pelstar® wireless technology |
| 2210CART | Optional rolling cart for 2210 series |
| SS-2210KL | ScaleSurance® 2 year extended warranty for 2210 series* |
| PELSTARONE-P | Wireless thermal printer |



2210CART
Optional Rolling Cart

To learn more, call **1-800-815-6615**
or visit **homscales.com/2210KL-AM-C**

Health o meter®
Professional

*Not available in all countries.