


# THE AUTOMATED CALIBRATOR

The Health o meter® Professional Scales automated calibrator, is a robust innovative calibration system that allows biomedical engineers and other medical technicians to quickly and safely calibrate medical patient scales without having to transport and stack physical calibration weights on the scale. Eliminating the cumbersome physical calibration weights from the calibration process reduces the time needed for routine scale maintenance and eliminates the safety risks associated with handling heavy dead weights. Traceable to the NIST, the automated calibration system records all calibration activity eliminating manual record keeping and simplifying audit activities.



**CALIBOT-AC**  
Acute Care Calibrator

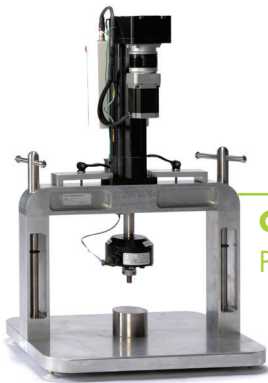

**CALIBOT-PC**  
Primary Care Calibrator

## ORDERING INFORMATION

MODEL	DESCRIPTION
CALIBOT-AC	Acute Care Calibrator
CALIBOT-PC	Primary Care Calibrator

Not available in all countries.

**Health o meter®**  
Professional

Health o meter® Professional Automated Calibrator

## TECHNICAL SPECIFICATIONS

### THE CALIBOT-AC ACUTE CARE CALIBRATOR

**Capacity**  
500 lb / 220 kg

**Resolution**  
0.01 lb / 0.005 kg

**Platform Size (w x d)**  
32" x 27"  
(813 mm x 656 mm)

**Product Weight**  
97 lb / 44 kg


### THE CALIBOT-PC PRIMARY CARE CALIBRATOR

**Capacity**  
500 lb / 220 kg

**Resolution**  
0.01 lb / 0.005 kg

**Platform Size (w x d)**  
17" x 17"  
(432 mm x 432 mm)

**Product Weight**  
65 lb / 30 kg

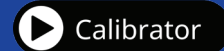

weigh**ea**sier®