

THE BCS-G7

BODY COMPOSITION SCALE

The BCS-G7 body composition scale features an industry-leading 1000 lb capacity, ensuring accurate weight and body composition measurements for any patient. With minimal pretest requirements, analysis can be performed at any time of day and under almost any circumstances, resulting in the greatest flexibility of any body composition analyzer. The large, spacious platform allows patients to find a comfortable and secure stance. Designed for patients 10 and over, the BCS-G7 delivers comprehensive body composition analysis, including detailed data for the torso and limbs, and differentiates between muscle mass and fat. This helps physicians track patient progress and response to treatment, especially for those using GLP-1 weight loss medications, ensuring patients are losing fat and increasing muscle mass.

The scale features a user-friendly interface making it the ideal body composition product for weight loss clinics, bariatric clinics, medspas, and doctor's offices that want easy functionality and fast results. The easy-to-read six-window display delivers rapid access to key body composition metrics including weight, body fat %, body water %, lean mass, Basal Metabolic Rate, and Body Mass Index.

The **BodyCompPro-G7** software simplifies weight tracking and enhances patient care with powerful, integrated features. Designed to connect directly to a computer, the scale's software enables seamless data tracking, detailed reporting, and streamlined patient management. Medical professionals benefit from easy-to-read summaries and graphic printouts—ideal for patient consultations and long-term monitoring.

THERE'S STRENGTH IN OUR NUMBERS.


When you add up the numbers on our scale, it's not surprising that the BCS-G7 body composition scale is a reliable, comprehensive choice for body composition analysis.

30 SECOND
TESTING TIME

11 BODY COMP
METRICS

120+ LOCAL
FIELD SPECIALISTS

100 YEARS
OF SERVICE

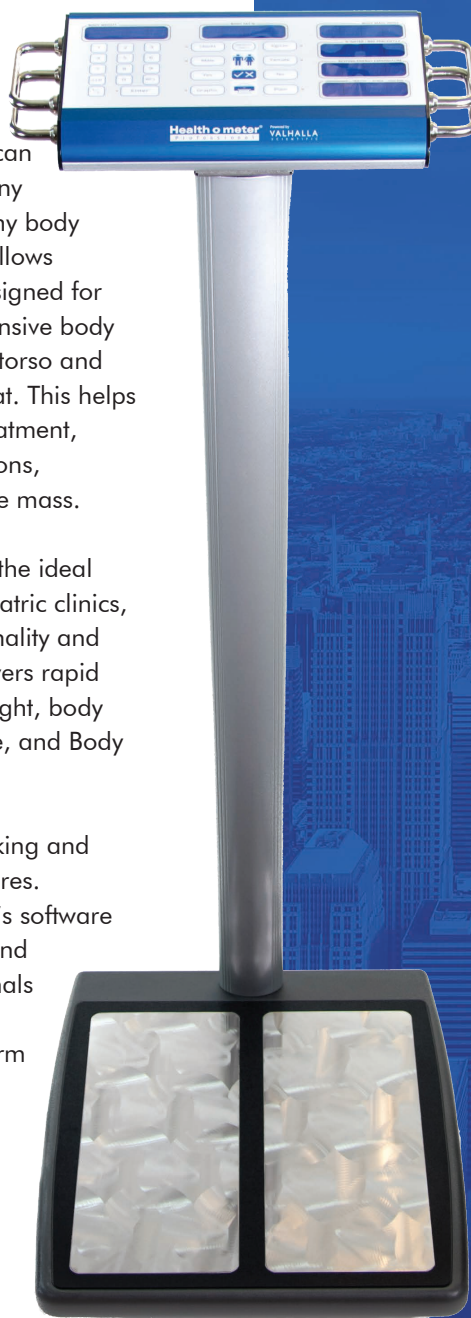

Health o meter®
Professional

Health o meter® Professional BCS-G7

weigh**easier**®

TECHNICAL SPECIFICATIONS

Capacity

1000 lb / 454 kg

Resolution

0.1 lb / 0.1 kg

Display

6-window display

Digit Size: 1/2"

Power Source

110-220V power adapter included

Connectivity

Integrated Software – Connect the scale to a computer for data tracking, reporting, patient management.*

Interface Type

Serial Interface to USB Converter

Functions

30-second testing start to finish

On Display Measurements:

Weight
Body fat %
Body water %
Muscle mass %
BMI
REE

Full Report Provides:

Weight
Body fat %
Body water %
Muscle mass %
Fat Free Mass
Ideal Weight Range
Segmental Fat Distribution Analysis
Total Body Fat
Total Body Water
BMI – Body Mass Index
BMR – Basal Metabolic Rate
REE – Resting Energy Expenditure
Active Metabolic Rate
Visceral Fat

Platform Size (w x d x h)

14" x 16" x 2 7/8"
(356 mm x 406 mm x 48 mm)

Product Footprint

20" x 18" x 41"
(508 mm x 457 mm x 1041 mm)

Product Weight

28 lb / 12 kg

Warranty

1-year limited manufacturer warranty

Bio-Impedance Specifications

Electrode Method: Tetrapolar System
Frequency: 50kHz
Test Current: 500µA
Conversion Rate: 750msec

ORDERING INFORMATION

MODEL	DESCRIPTION
BCS-G7	Body Composition Scale
IPO-ADU	Patient Report Template, Adult (5-page report, 25 reports/package)
IPO-ADU-SP	Spanish Language Patient Report Template, Adult (5-page report, 25 reports/package)
BCS-BROCHURE	Waiting Room Patient Brochure
SS-BCS-67	ScaleSurance® 2 year extended warranty for BCS-G7**

Scale includes integrated software. Visit homscscales.com/BCS-G7 to easily download BodyCompPro–G7: Body Composition App

To learn more, call **1-800-815-6615**
or visit **homscscales.com/BCS-G7**



*PC Requirements: Windows® 10/11 Operating System

**Not available in all countries.



Health o meter®
P r o f e s s i o n a l