

# The Weigh**easier**<sup>®</sup> Scale Selection Guide and Planning Tool

Designed for:

- Equipment Planners
- Architects
- Project Managers

**Health o meter**<sup>®</sup>  
P r o f e s s i o n a l



# why choose us?

**Pelstar LLC** is a leading manufacturer in the healthcare industry, known for its commitment to quality, innovation, and customer-focused solutions. As an ISO 13485 and 9001 certified company, Pelstar ensures its products adhere to the highest standards of safety and reliability. The company operates two renowned brands: **Health o meter<sup>®</sup> Professional Scales** and **Bridge Healthcare**.

What sets us apart is that everything we do - from product design to technical support - is done with our customers in mind. Each decision and each response is formed by first asking **What will be *weigh*easier for our customer?** – whether the customer is a physician, healthcare system, distributor, or sales team.

*weigh***easier<sup>®</sup>**

## #1 Medical Scale Brand in North America

Health o meter<sup>®</sup> has been a trusted brand since 1919. Today, Health o meter<sup>®</sup> Professional Scales is the U.S. market leader and the fastest growing brand in Medical Scales.\* We offer:

- High quality products that serve the entire continuum of medical care
- The expertise of MedTech/MedCare, the industry's largest network of trained and trusted field representatives
- Direct toll-free number for Technical Support with customer advocates based in the USA
- Fast shipments and industry-leading fill rates
- Standardization programs
- No-cost product evaluations

## 100% Focus on Healthcare

Because all we do is support the healthcare market, we understand the impact of healthcare reform, Value-Based Purchasing, and other industry-specific initiatives. We understand that providers are overwhelmed with patient care and are focused on operating efficiently and with an eye on costs. We concern ourselves with patient safety, supply chain efficiency, and technology standardization because our customers are concerned about those things. We design our products to help patients and the providers and institutions that care for them.

## Products Designed in the USA

Our US based engineers and product specialists design high-value, high-quality devices for users around the globe.

## Connectivity Solutions

Health o meter<sup>®</sup> Professional Scales provide seamless connectivity solutions for transmitting weight, height, and Body Mass Index (BMI) measurements to vital signs devices, PCs and EMR systems. Device connectivity enhances healthcare facility workflow efficiency, patient safety, and EMR compliance. Our compatibility extends to various systems, including:

WelchAlllyn<sup>®</sup>

SCHILLER

The Art of Diagnostics

MIDMARK<sup>®</sup>

EDANUSA

MDPRO

Allscripts<sup>®</sup>

Developer Program

CEIBA HEALTH

IOMT + VIRTUAL CLINICAL CARE

charmhealth

# Simplify Scale Selection with Confidence

The purpose of this Health o meter<sup>®</sup> Professional's Scale Selection Guide is to help ensure the right scales are chosen for a facility. We recognize that sorting through the options can be overwhelming, so we've created this guide to help narrow down the choices.

We also feature an interactive "Choose the Right Scale" tool on our website at: [www.homscales.com/products/choose-the-right-scale](http://www.homscales.com/products/choose-the-right-scale)

## Table of Contents

<div>Page 3</div> <div></div> <div>Scale Education</div>	<div>Page 4</div> <div></div> <div>ADA Compliant Scales</div>	<div>Page 5</div> <div></div> <div>Questions for Clients</div>	<div>Page 7</div> <div></div> <div>Acute Care Scales</div>
<div>Page 9</div> <div></div> <div>Primary Care Scales</div>	<div>Page 11</div> <div></div> <div>Most Popular Scales</div>	<div>Page 13</div> <div></div> <div>Stand-on Scales</div>	<div>Page 22</div> <div></div> <div>Wheelchair Scales</div>
<div>Page 27</div> <div></div> <div>Chair Scales</div>	<div>Page 28</div> <div></div> <div>Pediatric Scales and Carts</div>	<div>Page 32</div> <div></div> <div>Specialty Scales</div>	<div>Page 34</div> <div></div> <div>Height Measuring</div>
<div>Page 35</div> <div></div> <div>Scale Service and Maintenance</div>	<div>Page 36</div> <div></div> <div>Customer Support</div>		

1 \*As reported by Clarivate's US Medical Supply Distribution data.



# Scale Education

## Types of Scales

Health o meter® Professional Scales offers a variety of professional scales designed specifically to meet the needs of healthcare facilities.

### Stand-On Scales:

- Used by Primary Care Physicians, Clinics, and Hospitals
- Feature a variety of capacities and sizes, some with handrails
- Simple-to-use functions that make it user friendly and increases workflow efficiency
- Ideal for Pediatric, Adolescent, and Adult patients

### Wheelchair Scales:

- Used by Hospitals, Clinics, and Long Term Care facilities
- Feature high capacities, wide platforms, wheelchair ramps, and wheels for mobility
- Intuitive interface and functionality – including Tare which allows user to subtract the weight of the wheelchair
- Ideal for elderly, immobile, and bariatric patients

### Pediatric Scales:

- Used by Physician offices, Clinics, and Hospitals
- Feature easy-to-clean surfaces, built-in measuring tapes
- Effortless operations that allow for fast and easy weighing procedures

### Specialty Scales:

- Used by Hospitals, Clinics, and Long Term Care facilities
- Portfolio includes Diaper/Lap Sponge scales for OR, Transfer Board scale for ER/Stroke Units, Body Composition scales for Weight Loss Clinics

## Accessories/Options

After choosing the type of scale needed, there are variations and accessories to be considered.

- Handrails – Select scale models have integrated live handrails and some models have optional handrails to be added at the time of purchase or in the future.
- Remote Display – Scales with a remote display provide flexibility of allowing the display head to be mounted on a wall or placed on a table.
- Connectivity – Most digital scales have the option to connect to vital signs devices, PCs, and EMR systems via a variety of wireless kits. Selection of the appropriate kit requires a discussion between sales representative and the organization's IT team.
- Carts – Mobile carts are available to transport pediatric scales or other medical devices and supplies.
- KG Only – For facilities with metric requirements, scales can be ordered as KG models that provide weight measurements in kilograms only.
- ADA Compliant – Certified ADA compliant stand-on and wheelchair scales are available to help facilities meet U.S. Access Board requirements. See page 4 for more info.
- Height Rods – Stand-on scales may feature a built-in height rod or have the compatibility for an add-on height rod option. Wall-mounted models are also available.

# ADA Compliant Scales

At Health o meter® Professional Scales, safety and accessibility are at the forefront of our product design. We prioritize ADA compliance to ensure our healthcare equipment meets the highest standards. As the leading medical scale manufacturer in the U.S., we are dedicated to providing guidance and support as market requirements evolve.

In July 2024, the U.S. Access Board published new standards for medical diagnostic equipment (MDE) to ensure accessibility for all patients. Healthcare providers must comply with these standards within specified timelines:

- **Starting October 8, 2024:** All medical device equipment acquired by state and local government entities must be accessible.
- **By August 9, 2026:** Facilities must have at least 10% of their equipment, but no fewer than one unit, meeting the standards (20% for mobility rehabilitation facilities).

### Wheelchair Scale MDE Requirements:

- 1 32" wide by 40" deep platform
- 2 2" edge protection on platform
- 3 2" edge protection on ramps
- 4 1:8 running ramp slope



2101KL Digital Platform Scale



1100KL Digital Platform Scale with Extra-wide Handrails



### Stand-on Scale MDE Requirements:

- 1 Slip resistant platform
- 2 Handrail supports on both sides of platform
- 3 34" – 38" max. handrail height from platform

Learn More About our ADA Compliant Scales and the New MDE Standards



Source: U.S. Access Board. (2024). Medical Diagnostic Equipment (MDE) Accessibility Standards: <https://www.access-board.gov/mde/chapter-m3/>

# Questions for Clients

Health o meter® Professional's comprehensive line of scales addresses any need an organization might have with several options to consider. Here are some questions to discuss with an organization to help streamline the choices.

## Space Requirements

### ***Does the location of the scale placement have footprint limitations?***

- Example: If the specified location has a dimensional area of 36" x 50" then buyers should look closely at a scale's footprint size to ensure that the scale chosen will fit in the specified location, including surrounding area to allow for sufficient space to approach the scale. Models with foldable platforms are also available to allow for more flexible use of space.

## Portability

### ***Will the scale be moved throughout a department or stay in a central triage area?***

- Example: If the facility plans to move the scale frequently then models with portability features should be considered. Some models include transport handles and heavy-duty wheels or optional carts and travel bags.

## Platform Size and Capacity

### ***Will the patient population using the scale require a large platform and high capacity?***

- Example: A Bariatric Physician's office will require a scale with a larger platform to provide ample space for bariatric patients to have a comfortable and secure stance along with a high capacity to accommodate any size patient. Primary care clinics should also consider using large platform scales for elderly and unsteady patients, as these scales provide greater stability and are easier for them to stand on.

## Handrail Requirements

### ***Will the patient population using the scale require additional support of integrated handrails?***

- Example: Clinics with a majority of elderly patients would benefit from a scale that includes handrails. Equipment with "live" handrails allow patients to hold on for security giving the patient a more stable weighing experience.

## Unit Lock

### ***Are there any initiatives or requirements at the facility concerning measurement units?***

- Example: Some facilities require weight measurements to be recorded in kilograms. Most scales can be purchased as a KG only model so there is no risk of incorrect units being documented.

# Questions for Clients

## Workflow and Patient Population

### ***Does the facility require a scale that can accommodate both ambulatory patients and wheelchair users?***

- Example: If the facility sees a variety of patients, it would be beneficial to utilize a wheelchair scale. Ambulatory patients can stand on the scale platform during weighing and the same scale can be used to weigh patients in a wheelchair.

## ADA Requirements

### ***Has the facility considered the new ADA medical device requirements in their planning?***

- Example: By August 9, 2026 facilities must have at least 10% of their medical device equipment be accessible. These new requirements affect both equipment purchasing as well as planning for the space surrounding the equipment.

To learn more visit: <https://www.homscales.com/ada-mde/>

## EMR Connectivity

### ***Does the scale need the ability to connect to an EMR/EHR system?***


- Example: Some healthcare organizations have initiatives that require its medical devices to send patient data automatically into the facility's EMR system. These connectivity prerequisites require detailed conversations between the customer and the manufacturer and should be addressed when considering scale purchases.

## Height Measuring

### ***Will the provider need a way to measure patient height?***

- Example: Primary care clinics often measure patient height in conjunction with weight to determine Body Mass Index or evaluate pediatric patient's growth.

# Acute Care Scales

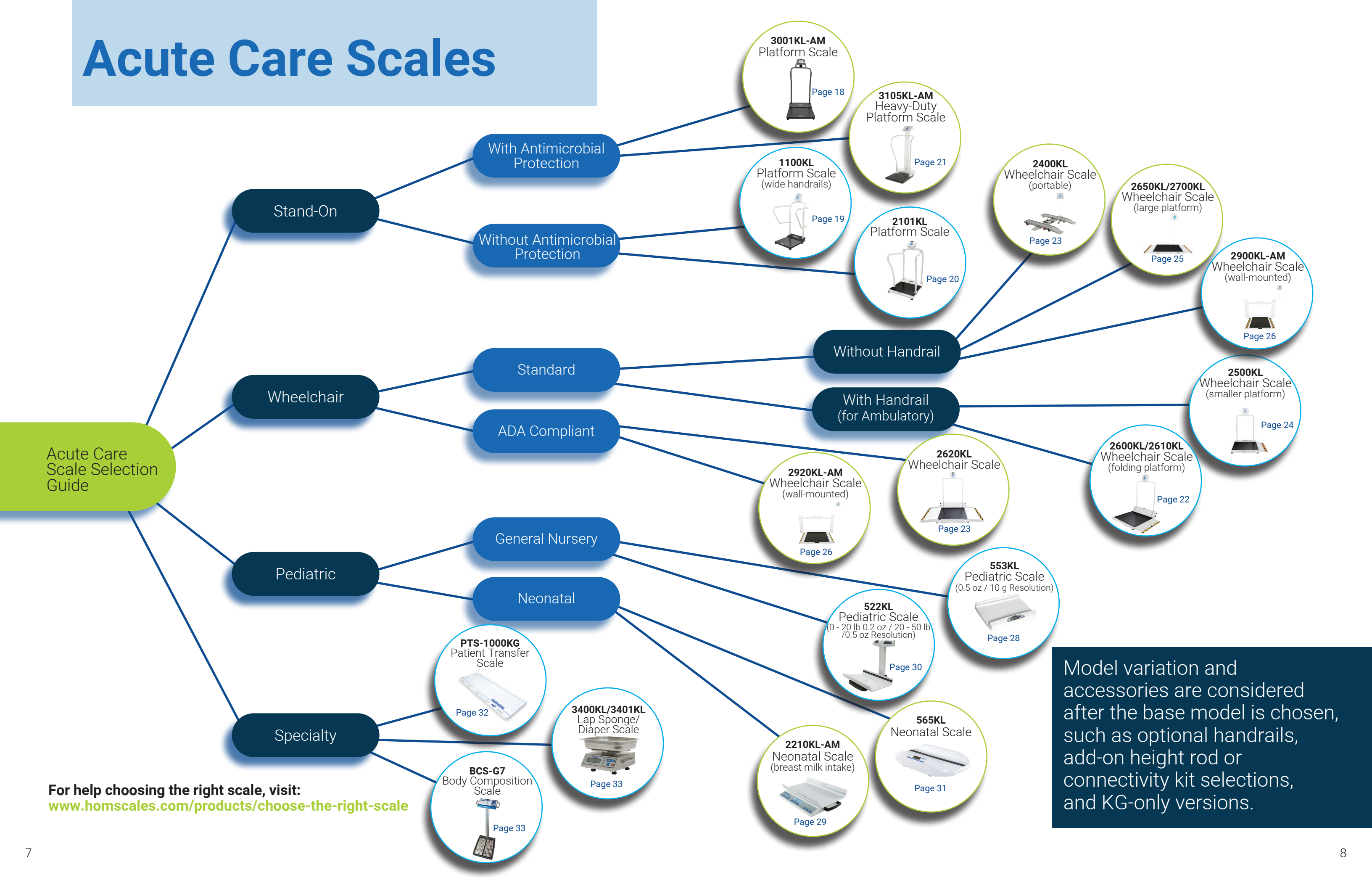


For help choosing the right scale, visit:  
[www.homscales.com/products/choose-the-right-scale](http://www.homscales.com/products/choose-the-right-scale)

Model variation and accessories are considered after the base model is chosen, such as optional handrails, add-on height rod or connectivity kit selections, and KG-only versions.



# Primary Care Scales

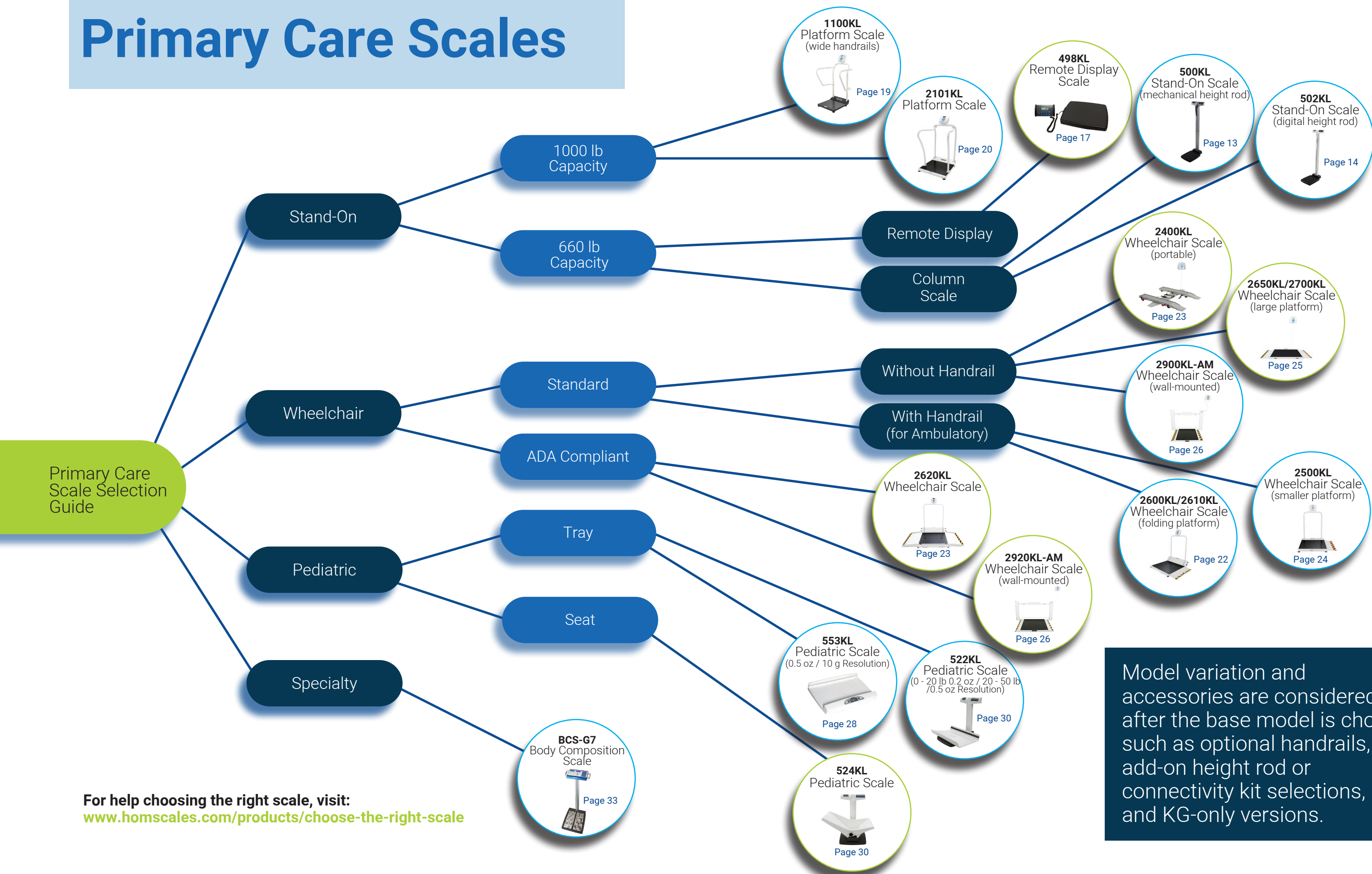

For help choosing the right scale, visit:  
[www.homscales.com/products/choose-the-right-scale](http://www.homscales.com/products/choose-the-right-scale)

Model variation and accessories are considered after the base model is chosen, such as optional handrails, add-on height rod or connectivity kit selections, and KG-only versions.

# Most Popular Acute Care Scales



## Heavy Duty Mobile Antimicrobial Platform Scale

3105KL-AM

- Capacity of 1000 lb/ 454 kg
- EMR connectivity, BMI, hold/release, reweigh, recall, zero, KG/LB, tare, pre-tare, Everlock®
- Smooth-rolling wheels and large handle make transport easy
- Optional electronic and mechanical height rod
- Improved robust structural design and low profile platform
- Options such as KG only and built-in wireless models are available



## High-Resolution Antimicrobial Digital Neonatal/Pediatric Tray Scale

2210KL-AM (2-sided tray)  
2210KL4-AM (4-sided tray)  
2210CART-AC (acute care cart)

- High 1 gram resolution with a capacity of 20,000g / 20kg / 45lb
- EMR connectivity, breast milk intake, auto-weight hold, hold/release, zero, tare, pre-tare, reweigh, recall, GM/ KG/LB conversion and lockout, auto-zero, and auto-off
- Round-edged removable tray and backlit display
- Optional electronic height rod
- Options such as KG only, GM only, LB only, and built-in wireless models are available




## Wall-Mounted Fold-Up Antimicrobial Wheelchair Scale

2900KL-AM

2920KL-AM (ADA Compliant)

- Capacity of 1000 lb/454 kg
- Large easy to lift platform conserves space when folded
- EMR connectivity, BMI, hold/release, reweigh, recall, zero, KG/ LB, tare, pre-tare, Everlock®
- Safety features include slip resistant caution tape, recessed handle, raised edges, and no pinch points
- Options such as KG only, and built-in wireless models are available
- ADA Compliant model available to support facilities in meeting accessibility standards

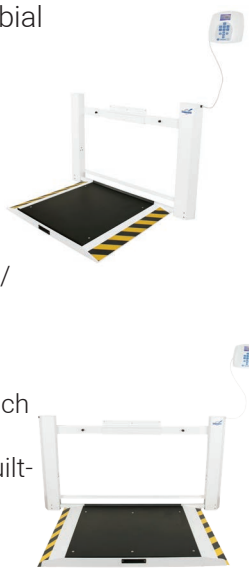

## Patient Transfer Scale

PTS-1000KG

- A transfer board with a built-in medical scale that provides a weight measurement when transferring a patient
- Weighs immobile patients instantly and accurately
- It helps shorten door-to-needle times for time-critical patients
- All departments can benefit from the PTS, particularly Stroke Units, the Emergency Room, ICU, and Radiology
- KG only model is a safe solution for facilities standardizing on metric only measurements for critical medication dosages based on kilograms



## Digital Wet Diaper/Lap Sponge/ Organ Scale

3400KL / 3400KG (Pan)

3401KL (Tray)

- Stable, stainless steel base and removable stainless steel weighing pan or tray that can be autoclaved for sterilization
- Powered by a long-life rechargeable battery
- Backlit display that's easy to read at any angle and light condition



# Most Popular Primary Care Scales



## Primary Care Digital Scale

500KL/ 502KL (digital height rod)

- Capacity of 660 lb/300 kg
- BMI, hold/release, reweigh, recall, zero, KG/LB, tare, pre-tare, Everlock® and more
- 8" wide wheelbase for easy portability
- The most extensive connectivity options in the market
- Options such as KG only and built-in handlebars



## Remote Display Digital Scale

498KL

- Capacity of 660 lb/300 kg
- BMI, hold/release, reweigh, recall, zero, KG/LB, tare, pre-tare, Everlock® and more
- Remote display for easy portability
- Large non-skid platform with minimal step-up
- Options such as KG only



## Digital Pediatric Scale with Extra-Wide Tray

553KL

- Capacity of 44 lb/20 kg
- Hold/release, reweigh, recall, zero, KG/ LB, tare, pre-tare, Everlock® and more
- Measuring tape included
- Extra-Wide tray provides additional patient comfort and safety
- Optional electronic and mechanical height rod
- Motion Detection – accurately weigh patients even in the presence of movement
- Options such as KG only



## Digital Pediatric Tray Scale with Raised Display

522KL

- Capacity of 44 lb/20 kg
- Hold/release, reweigh, recall, zero, KG/ LB, tare, pre-tare, Everlock® and more
- Measuring tape included
- Elevated display head for easy-to-read weight measurements
- Optional electronic and mechanical height rod
- Motion Detection – accurately weigh patients even in the presence of movement
- Options such as KG only





# Stand-On Scales



★TOP SELLER  
**500KL**



Videos & More!

EYE-LEVEL DIGITAL SCALE

The 500KL is the #1 selling physician scale in North America. With a built-in height rod, the scale enhances workflow by measuring weight and height and calculating Body Mass Index without moving the patient. Optional handlebars provide stability to keep both patients and caregivers safe.

**Built-in Height Rod:** The easy-to-use height rod provides height measurements from 30" - 84 1/4" / 76 cm - 214 cm.

**Calculates BMI:** Scale functions include the ability to calculate the patient's BMI.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.

**Dent/Rust Proof Low Profile Platform:** Non-skid platform requires only minimal step up, offering the patient a comfortable and stable weighing experience.



Capacity	Resolution	Display	Platform Size (w x d x h)
660 lb / 300 kg	0.2 lb / 0.1 kg	1" / 25 mm LCD	13 3/4" x 16 1/2" x 2 3/8" 349 mm x 419 mm x 60 mm

Height Rod	Height Rod Graduation	Product Footprint (w x d x h)
30" - 84 1/4" 76 cm - 214 cm	1/8" / 1 mm	14 3/8" x 21 1/8" x 53 5/8" 359 mm x 537 mm x 1362 mm

#### Options & Accessories:

Model	Description
500KG	Scale in kilograms only
500KLHB	500KL with handlebars included
500KLAD	500KL with adapter included
500KGAD	500KG with adapter included
500KL-BT	500KL with built-in Pelstar® wireless technology
500KG-BT	500KG with built-in Pelstar® wireless technology
500HB	Optional handlebar accessory for 500 series scales
SS-500KL	2 year extended warranty
ADPT31	Adapter for 500 series scales
BTPRINTERKIT-01	Wireless thermal printer kit

**Power Source:** 6 AA batteries included, 100-240V adapter (not included)

**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Auto Zero, Auto Off, Reweigh, Recall, Parent / Child Weighing Function, Variable Auto Off Time, Audible / Mute Sound Option, Enabling Unit View with Unit of Measure Locked, Everlock®

**Wheels:** 2 wheels, 8" wide wheelbase

**Connectivity:** optional Pelstar® wireless technology



# 502KL



Videos & More!

PRIMARY CARE DIGITAL SCALE  
WITH DIGITAL HEIGHT ROD

The newest member of the #1 selling physician scale in North America. With a market-leading built-in digital height rod, the scale enhances workflow by measuring weight and height and calculating Body Mass Index without moving the patient. Optional handlebars provide stability to keep both patients and caregivers safe.

**Built-in Digital Height Rod:** The easy-to-use height rod provides a full measuring range of 11 3/4" - 79 1/2" / 30cm - 202 cm and a digital range of 3ft 11.2" - 6ft 7.5" / 120 cm - 202 cm.

**Calculates BMI:** Scale functions include the ability to calculate the patient's BMI.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.

**Dent/Rust Proof Low Profile Platform:** Non-skid platform requires only minimal step up, offering the patient a comfortable and stable weighing experience.



Optional 500HB  
Handlebar Accessory  
can be mounted to  
your 502KL,  
502KG,  
502KL-BT, and  
502KG-BT



Capacity	Resolution	Display	Platform Size (w x d x h)
660 lb / 300 kg	0.2 lb / 0.1 kg	1" / 25 mm LCD	13 3/4" x 16 1/2" x 2 3/8" 349 mm x 419 mm x 60 mm

Digital Height Rod	Mechanical Height Rod	Product Footprint (w x d x h)
3ft 11.2" - 6ft 7.5" 120 cm - 202 cm	11 3/4" - 79 1/2" 30 cm - 202 cm	14 3/8" x 21 1/8" x 53 5/8" 359 mm x 537 mm x 1362 mm
Resolution: 1 mm / 0.1"		Graduation: 1/16" / 1 mm

#### Options & Accessories:

Model	Description
502KG	Scale in kilograms only
502KL-BT	502KL with built-in Pelstar® wireless technology
502KG-BT	502KG with built-in Pelstar® wireless technology
500HB	Optional handlebar accessory for 502 series scales
SS-502KL	2 year extended warranty
ADPT31	Adapter for 502 series scales
BTPRINTERKIT-01	Wireless thermal printer kit

**Power Source:** 6 AA batteries included, 100-240V adapter (not included)

**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Auto Zero, Auto Off, Reweigh, Recall, Parent / Child Weighing Function, Variable Auto Off Time, Audible / Mute Sound Option, Enabling Unit View with Unit of Measure Locked, Everlock®

**Wheels:** 2 wheels, 8" wide wheelbase

**Connectivity:** optional Pelstar® wireless technology



499KL

WAIST-HIGH DIGITAL SCALE



The 499KL has a waist-high easy to read display and the scale is an attractive, robust option for use in medical offices. With an optional height rod and optional handlebars, the 499KL is an economical choice for almost any environment.

**Capacity:** 460 lb / 209 kg

**Resolution:** 0.2 lb / 0.1 kg

**Platform Size** (w x d x h):

13 3/4" x 16 1/2" x 2 3/8"

349 mm x 419 mm x 60 mm

**Wheels:** 2 wheels, 8" wide wheelbase

**Product Footprint** (w x d x h):

14 1/8" x 21 1/8" x 43 5/8"

359 mm x 537 mm x 1108 mm

**Options & Accessories:**

**Display:** 1" / 25 mm LCD

**Power Source:** 6 AA batteries

included, 100-240V adapter

(not included)

**Functions:** LB / KG Conversion,  
LB / KG Lock Out, BMI, Zero, Tare,  
Hold / Release, Auto Zero, Auto Off,  
Parent / Child Weighing Function,  
Everlock®

Model	Description
499KG	Scale in kilograms only
499KLHR	499KL with height rod included
499KGHR	499KG with height rod included
499KLAD	499KL with adapter included
499KGAD	499KG with adapter included
499KLHRAD	499KL with height rod and adapter included
499KGHRAD	499KG with height rod and adapter included
499HB	Optional handlebar accessory for 499 series scales
SS-499KL	2 year extended warranty
499KLROD	Height rod for 499 series scales
ADPT31	Adapter for 499 series scales



599KL

HEAVY DUTY WAIST-HIGH DIGITAL SCALE



Great for doctor offices and clinics with minimal floor space, the 599KL has a high capacity with waist-high display that can swivel 350°. For providers wishing to measure height, an optional height rod is available.

**Capacity:** 660 lb / 300 kg

**Resolution:** 0.2 lb / 0.1 kg

**Platform Size** (w x d x h):

14 1/4" x 14 1/4" x 2 5/8"

362 mm x 362 mm x 67 mm

**Connectivity:** USB, optional Pelstar  
wireless technology

**Product Footprint** (w x d x h):

14 3/4" x 23 1/4" x 42 1/2"

375 mm x 591 mm x 1384 mm

**Options & Accessories:**

**Display:** 1" / 25 mm LCD, 350°

Swivel

**Wheels:** 2

**Power Source:**

100-240V adapter included,

6 C batteries (not included)

**Functions:** LB / KG Conversion,  
LB / KG Lock Out, BMI, Zero, Tare,  
Hold / Release, Auto Zero, Auto Off,  
Parent / Child Weighing Function

Model	Description
599KG	Scale in kilograms only
599KLHR	599KL with height rod included
599KGHR	599KG with height rod included
SS-599KL	2 year extended warranty
499KLROD	Height rod for 599KL
ADPT50	Adapter included, item # for replacement or spare



600KL

HEAVY DUTY EYE-LEVEL DIGITAL SCALE  
WITH DIGITAL HEIGHT ROD

Improve workflow with this top-of-the-line physician office scale. The built-in digital height rod eliminates the need to manually enter height for computing BMI. The 350° swivel head allows for line-of-sight viewing and patient privacy.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.

**Automatically Calculates Weight, Height and BMI:** Improve efficiency and accuracy of vital sign measurements.

**Built-in Digital Height Rod:** The easy-to-use, state of the art digital height rod provides measurements from 23 5/8" - 88" / 60 cm - 222 cm.

**Patented Design:** The 600KL's patented design provides a user-friendly and intuitive interface for measuring patient weight, height and BMI with a convenient rotating display head.



Capacity	Resolution	Platform Size (w x d x h)
----------	------------	---------------------------

660 lb / 300 kg    0.2 lb / 0.1 kg    14 1/4" x 14 1/4" x 2 5/8"  
362 mm x 362 mm x 67 mm

Height Rod	Display	Product Footprint (w x d x h)
------------	---------	-------------------------------

23 5/8" - 88"  
60 cm - 222 cm  
Resolution: 0.1" / 0.1 cm    1" / 25 mm  
LCD 350° Swivel    14 3/4" x 24 1/8" x 54 3/8"  
375 mm x 613 mm x 1381 mm

**Options & Accessories:**

Model	Description
600KG	Scale in kilograms only
600KL-BT	600KL with built-in Pelstar® wireless technology
600KG-BT	600KG with built-in Pelstar® wireless technology
ADPT50	Adapter included, item # for replacement or spare
SS-600KL	2 year extended warranty

498KL

REMOTE DISPLAY DIGITAL SCALE



The 498KL is a versatile choice for almost any clinical environment. Its portability and remote display head allow the scale to be moved from one location to another.

**Capacity:** 660 lb / 300 kg  
**Resolution:** 0.2 lb / 0.1 kg  
**Platform Size (w x d x h):**  
13 ¾" x 16 ½" x 2 ⅝"  
349 mm x 419 mm x 60 mm  
**Connectivity:** via optional Pelstar® wireless technology  
**Product Footprint (w x d x h):**  
13 ¾" x 16 ½" x 2 ⅝"  
349 mm x 419 mm x 60 mm

Options & Accessories:

Model	Description
498KG	Scale in kilograms only
498KLAD	498KL with adapter included
SS-498KL	2 year extended warranty
ADPT31	Adapter for 498KL
64771	Carrying case for 498KL
BTPRINTERKIT-01	Wireless thermal printer kit



**Display:** 1" / 25 mm LCD  
**Power Source:** 6 AA batteries included, 100-240V adapter (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Reweigh, Auto Zero, Auto Off, Parent / Child Weighing Function, Variable Auto Off Time, Audible / Mute Sound Option, Enabling Unit View with Unit of Measure Locked, Everlock®

752KL

HEAVY DUTY REMOTE DISPLAY DIGITAL SCALE



With a high 660 lb capacity, the 752KL is a great supplement to accommodate ambulatory bariatric patients. The remote display and portable base make it convenient to move the scale anywhere it is needed.

**Capacity:** 660 lb / 300 kg  
**Resolution:** 0.2 lb / 0.1 kg  
**Platform Size (w x d x h):**  
14 ¼" x 14 ¼" x 2 ⅝"  
362 mm x 362 mm x 67 mm  
**Connectivity:** USB, optional Pelstar® wireless technology  
**Product Footprint (w x d x h):**  
14 ½" x 14 ½" x 2 ¾"  
368 mm x 368 mm x 70 mm

Options & Accessories:

Model	Description
752KG	Scale in kilograms only
SS-752KL	2 year extended warranty
ADPT50	Adapter included, item # for replacement or spare
64771	Carrying case for 752KL



**Display:** 1" / 25 mm LCD  
**Power Source:** 100-240V adapter included, 6 C-cell batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Hold / Release, Auto Zero, Auto Off, Parent / Child Weighing Function

Optional Digital and Mechanical Height Rods Available



3001KL-AMX  
ANTIMICROBIAL DIGITAL  
PLATFORM SCALE WITH  
EXTENDED HANDRAILS

**Power Source:** 100-240V adapter included, 6 D batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Audible / Mute Sound Option, Parent / Child Weighing Function, Everlock®  
**Wheels:** 2  
**Connectivity:** USB, optional Pelstar® wireless technology

★TOP SELLER

3001KL-AM  
& 3001KL-AMX

ANTIMICROBIAL DIGITAL PLATFORM SCALE

Incredibly robust, these platform scales feature feet and wheels specially designed to withstand frequent use and regular movement. Antimicrobial coating on the handrails, transport handle, and keypad help reduce the spread of harmful microbes. The 3001KL-AMX offers extended handrails that provide extra support and help patients maintain balance when stepping onto the platform.

**Infection Prevention:** Silver ionic Antimicrobial Powder Coating suppresses the growth of bacteria and helps protect against cross-contamination.

**Fall Prevention:** "Live" handrails are advanced technology that allow the patient to hold on for stability while the scale accurately calculates weight.

**Low Profile Platform:** Requires only minimal step up, offering the patient a comfortable and stable weighing experience.

**Durability:** Convenient transport handle allows the user to easily move the scale to where it is needed in the facility. Rugged wheels allow for smooth and easy portability.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.



Capacity	Resolution	Platform Size (w x d x h)
1000 lb / 454 kg	0.2 lb / 0.1 kg	22" x 15 ¾" x 2 ⅝" 559 mm x 400 mm x 67 mm

Display	Product Footprint (w x d x h)
Digit Size: ¾" / 10 mm, Screen Size: 2 ¾" x 1 ½" High-Contrast Color TFT-LCD Screen	22" x 33 ⅞" X 52" 559 mm x 841 mm x 1321 mm

Options & Accessories:

Model	Description
3001KG-AM	Assembled scale in kilograms only
3001KL-AM-BT	3001KL-AM with built-in Pelstar® wireless technology
3001KG-AM-BT	3001KG-AM with built-in Pelstar® wireless technology
3001KL-AMHR	3001KL-AM with mechanical height rod
3001KG-AMHR	3001KG-AM with mechanical height rod
3001KL-AMX	Assembled scale with extended handrails
3001KG-AMX	3001KL-AMX in kilograms only
245EHR-3001	Digital height rod for 3001KL-AM & 3001KG-AM
3001HR	Mechanical height rod for 3001KL-AM & 3001KG-AM
SS-3001KL-AM	2 year extended warranty
SS-3001KL-AMX	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit





Optional Digital and Mechanical Height Rods Available



★ TOP SELLER

# 1100KL

## DIGITAL PLATFORM SCALE WITH EXTRA WIDE HANDRAILS

This Certified ADA Compliant platform scale supports up to 1000 pounds, and offers “live” handrails on both sides of the slip resistant low profile platform to provide a stable and comfortable weighing experience for patients of any size. The keypad features antimicrobial protection. Ideal for any type of facility including primary care, rehab or hospital.

**Fall Prevention:** “Live” handrails are advanced technology that allow the patient to hold on for stability while the scale accurately calculates weight.

**Low Profile Platform:** Requires only minimal step up, offering the patient a comfortable and stable weighing experience.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.

**Complies with MDE Standards:** Slip resistant platform, handrail supports on both sides of platform, and 34” min. – 38” max. handrail height from platform



Capacity	Resolution	Platform Size (w x d x h)
1000 lb / 454 kg	0.2 lb / 0.1 kg	22" x 15 ¾" x 2 ½" 559 mm x 400 mm x 64 mm

Display	Product Footprint (w x d x h)
---------	-------------------------------

Digit Size: ¾" / 10 mm, Screen Size: 2 ¾" x 1 ½" High-Contrast Color TFT-LCD Screen	31 ¾" x 33 ½" x 49 ¾" 806 mm x 851 mm x 1264 mm
---	--

### Options & Accessories:

Model	Description
1100KG	Scale in kilograms only
1100KL-BT	1100KL with built-in Pelstar® wireless technology
1100KG-BT	1100KG with built-in Pelstar® wireless technology
1100KLHR	1100KL with mechanical height rod included
1100KGHR	1100KG with mechanical height rod included
1100KL-EHR	1100KL with digital height rod included
1100KL-BT-EHR	1100KL-EHR with built-in Pelstar® wireless technology
245EHR-1100	Digital height rod for 1100 series scales
PROPLUSROD	Mechanical height rod for 1100 series scales
SS-1100KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit

**Power Source:** 100-240V adapter included, 6 D batteries (not included)

**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Audible / Mute Sound Option, Parent / Child Weighing Function, Everlock®

**Wheels:** 2  
**Connectivity:** USB, optional Pelstar® wireless technology

## 2101KL

### DIGITAL PLATFORM SCALE



Optional Digital and Mechanical Height Rods Available



This Certified ADA Compliant bariatric scale features “live” handrails on both sides of the platform that allow the patient to hold on for stability while taking an accurate weight. The keypad features antimicrobial protection. Suitable for use in MD offices, long term care, hospitals or anywhere there is adequate space for the device.

- Capacity:** 1,000 lb / 454 kg

**Resolution:** 0.2 lb / 0.1 kg

**Platform Size (w x d x h):**  
26" x 22" x 2 ¾"  
660 mm x 559 mm x 60 mm

**Display:** Digit Size: ¾" / 10 mm,  
Screen Size: 2 ¾" x 1 ½" High-Contrast  
Color TFT-LCD Screen

**Product Footprint (w x d x h):**  
31" x 27" X 50"  
787 mm x 686 mm x 1270 mm

**Options & Accessories:**
- Connectivity:** Via USB, optional Pelstar® wireless technology

**Power Source:** 100-240V adapter included, 6 D batteries (not included)

**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Audible / Mute Sound Option, Parent / Child Weighing Function, Everlock®

Model	Description
2101KG	Scale in kilograms only
2101KL-BT	2101KL with built-in Pelstar® wireless technology
2101KG-BT	2101KG with built-in Pelstar® wireless technology
2101KLHR	2101KL with mechanical height rod included
2101KGHR	2101KG with mechanical height rod included
245EHR-2101	Digital height rod for 2101 series scales
STROD	Mechanical height rod for 2101 series scales
SS-2101KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit



## 1110KL

### HIGH CAPACITY REMOTE DISPLAY DIGITAL SCALE WITH LARGE BASE



Optional Accessories Include Handlebars, and a Digital or Mechanical Height Rod

A versatile high-capacity remote digital scale designed to be customizable to accomodate any patient or facility. The keypad features antimicrobial protection.

- Capacity:** 1000 lb / 454 kg

**Resolution:** 0.2 lb / 0.1 kg

**Platform Size (w x d x h):**  
22" x 15 ¾" (d) x 2 ⅝"  
559 mm x 597 mm x 121 mm

**Display:** Digit Size: ¾" / 10 mm,  
Screen Size: 2 ¾" x 1 ½" High-Contrast  
Color TFT-LCD Screen

**Product Footprint (w x d x h):**  
22" x 23 ½" X 4 ¾"  
559 mm x 597 mm x 121 mm

**Options and Accessories:**
- Connectivity:** Via USB, optional Pelstar® wireless technology

**Power Source:** 100-240V adapter included, 6 D batteries (not included)

**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Audible / Mute Sound Option, Parent / Child Weighing Function, Everlock®

Model	Description
1110KG	Scale in kilograms only
1110KL-BT	1110KL with built-in Pelstar® wireless technology
1110KG-BT	1110KG with built-in Pelstar® wireless technology
1110HB	Handlebars for the 1110 series scales
245EHR-1110	Digital height rod for 1110 series scales (must have an 1110HB)
201HR-1110	Mechanical height rod for 1110 series scales (must have an 1110HB)
SS-1110KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit





Optional Digital and Mechanical Height Rods Available



**PowerSource:** 6 D batteries included, 100-240V adapter not included  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Audible / Mute Sound Option, Parent / Child Weighing Function, Everlock®  
**Connectivity:** USB, optional Pelstar® wireless technology

★TOP SELLER

# 3105KL-AM

HEAVY-DUTY ANTIMICROBIAL PLATFORM SCALE

The 3105KL-AM has a robust, innovative design with enhanced safety and mobility features, "live" handrails allow the patient to hold onto handrails for stability while the scale accurately calculates weight. Low profile platform requires only minimal step up, offering the patient a comfortable and stable weighing experience.

**Safety Features:** Low-profile platform, "live" handlebars, antimicrobial protection on the head, keypad, transport handle, handrails, body, and kickplate, and an innovative design that eliminates places for dirt and contaminants.

**Portability:** 2 smooth-rolling wheels and easy grip handle make it easy to move from room to room.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.



Capacity	Resolution	Platform Size (w x d x h)
1000 lb / 454 kg	0.2 lb / 0.1 kg	19 3⁄4" x 21 1⁄4" x 2 1⁄8" 502 mm x 540 mm x 54 mm
Display		Product Footprint (w x d x h)
Digit Size: 3⁄8" / 10 mm, Screen Size: 2 3⁄4" x 1 1⁄2" High-Contrast Color TFT-LCD Screen		23 1⁄4" x 28" x 52 1⁄4" 591 mm x 711 mm x 1327 mm

Options & Accessories:

Model	Description
3105KG-AM	Scale in kilograms only
3105KL-AM-BT	3105KL-AM with built-in Pelstar® wireless technology
3105KG-AM-BT	3105KG-AM with built-in Pelstar® wireless technology
245EHR-3105	Digital height rod for 3105KL-AM series scales
201HR-3105	Mechanical height rod for 3105KL-AM series scales
SS-3105KL-AM	2 year extended warranty
ADPT30	Adapter for 3105KL-AM series scales
BTPRINTERKIT-01	Wireless thermal printer kit

# Wheelchair Scales



Optional **2600ADAKIT** ADA Compliant Upgrade Kit can be installed on your existing 2600KL and 2610KL series of scales to bring your scale into compliance with the new MDE standards

**Power Source:** 100-240V adapter included, 6 D batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Everlock®, Audible / Mute Sound Option  
**Wheels:** 4, including 2 casters  
**Connectivity:** USB, optional Pelstar® wireless technology

★TOP SELLER

# 2600KL / 2610KL

DIGITAL WHEELCHAIR RAMP SCALE / DIGITAL WHEELCHAIR DUAL RAMP SCALE

Our most popular wheelchair scale is also the most versatile, and is a good choice for any facility. The platform folds for convenient storage and the casters (the only scale in the market that has them!) make transport a breeze. For easy on/off access the 2610KL offers dual ramps.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.

**Casters and Folding Ramps for Easy Portability and Storage:** 4 wheels, including 2 casters, afford ease of movement and the folding ramps make it simple to store.

**Wide Platform for Wheelchair or Bariatric Use:** Wide platform accommodates wheelchairs of all sizes and provides a comfortable weighing experience for bariatric patients.

**Calculates BMI:** Scale functions include the ability to calculate patient's BMI.



Capacity	Resolution	Platform Size (w x d x h)
1000 lb / 454 kg	0.2 lb / 0.1 kg	36" x 32 1⁄4" x 2 1⁄2" 914 mm x 819 mm x 64 mm
Display		Product Footprint (w x d x h)
Digit Size: 3⁄8" / 10 mm, Screen Size: 2 3⁄4" x 1 1⁄2" High-Contrast Color TFT-LCD Screen		48 1⁄8" x 43 7⁄8" x 49 3⁄8" (2600KL) 1222 mm x 1114 mm x 1254 mm 56 3⁄4" x 43 7⁄8" x 49 3⁄8" (2610KL) 1441 mm x 1114 mm x 1254 mm

Options & Accessories:

Model	Description
2600KG	Scale in kilograms only
2610KG	Scale in kilograms only
2600KL-BT	2600KL with built-in Pelstar® wireless technology
2610KL-BT	2610KL with built-in Pelstar® wireless technology
2600KG-BT	2600KG with built-in Pelstar® wireless technology
2610KG-BT	2610KG with built-in Pelstar® wireless technology
2620KL	ADA compliant scale
2620KG	ADA compliant scale in kilograms only
SS-2600KL	2 year extended warranty
SS-2610KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
2600ADAKIT	ADA Compliance Kit for gray 2600KL & 2610KL series scales
2600ADAKIT-WHITE	ADA Compliance Kit for white 2600KL & 2610KL series scales
B2600RAMP	Secondary ramp in white for 2600KL & 2600KG series scales
BTPRINTERKIT-01	Wireless thermal printer kit



2620KL

DIGITAL WHEELCHAIR SCALE WITH DUAL RAMP AND X-LARGE PLATFORM, ADA COMPLIANT



**Power Source:** 100-240V adapter included, 6 D batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Everlock®, Audible / Mute Sound Option  
**Wheels:** 4, including 2 casters  
**Connectivity:** USB, optional Pelstar® wireless technology

2400KL

DIGITAL PORTABLE WHEELCHAIR SCALE



Designed to meet the MDE Standards for Accessible Medical Diagnostic Equipment and is 3rd party independently certified in compliance with accessibility codes and standards. The scale includes critical features that meet ADA safety requirements and enhanced safety features to ensure a comfortable weighing experience for patients who use wheelchairs.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.

**Complies with MDE Standards:** Dual ramps adhere to a 1:8 slope ratio. The ramps and extra-large 32" x 40" platform feature a 2" edge protection.

**Safety Features:** Ramps and extra-large platform feature edge protection to prevent roll-offs, and anti-slip caution tape on the ramps enhances safety providing both the clinician and patient confidence during the weighing experience.



Capacity	Resolution	Platform Size
1000 lb / 454 kg	0.2 lb / 0.1 kg	40" x 32 ¼" x 1 ⅞" 1016 mm x 819 mm x 48 mm
Display		Product Footprint (w x d x h)
Digit Size: ⅜" / 10 mm, Screen Size: 2 ¾" x 1 ½" High-Contrast Color TFT-LCD Screen		77 ¼" x 43 ⅞" x 48 ¾" 1962 mm x 1114 mm x 1238 mm

Options & Accessories:

Model	Description
2620KG	Scale in kilograms only
SS-2620KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit

The 2400KL is the perfect scale for locations where occasional wheelchair-bound patients need to be weighed. The scale folds easily and requires minimal room for storage. The scale includes a handle for portability.

**Capacity:** 800 lb / 363 kg  
**Resolution:** 0.2 lb / 0.1 kg  
**Platform Size:** Expands and adjusts from 25 ¾" - 50 ⅞" / 654 mm - 1273 mm  
**Display:** Digit Size: ⅜" / 10 mm, Screen Size: 2 ¾" x 1 ½" High-Contrast Color TFT-LCD Screen  
**Product Footprint (w x d x h):** 51 ¼" x 40 ½" x 3 ½" (ramps fully extended)  
1302 mm x 1029 mm x 89 mm

Options & Accessories:

Model	Description
2400KG	Scale in kilograms only
2400KL-BT	2400KL with built-in Pelstar® wireless technology
2400KG-BT	2400KG with built-in Pelstar® wireless technology
SS-2400KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit

**Connectivity:** USB, optional Pelstar® wireless technology  
**Wheels:** 2 wheels, 2 casters  
**Power Source:** 100-240V adapter included, 6 D batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Everlock®, Audible / Mute Sound Option



2500KL

DIGITAL WHEELCHAIR RAMP SCALE



2500CKL

DIGITAL WHEELCHAIR RAMP SCALE WITH FOLD AWAY SEAT



This single-ramp wheelchair scale is perfect for use with standard size wheelchairs. To weigh patients in a seated position, a chair accessory is also available, item # CHAIRACC.

**Capacity:** 1,000 lb / 454 kg  
**Resolution:** 0.2 lb / 0.1 kg  
**Platform Size:** 29 ⅜" (w) x 28 ¼" (d) x 2 ⅝" (h)  
746 mm (w) x 718 mm (d) x 60 mm (h)  
**Display:** Digit Size: ⅜" / 10 mm  
Screen Size: 2 ¾" x 1 ½" High-Contrast Color TFT-LCD Screen  
**Connectivity:** USB, optional Pelstar® wireless technology  
**Product Footprint (w x d x h):** 40 ¼" x 34 ½" x 45 ½"  
1022 mm x 876 mm x 1156 mm  
**Options and Accessories:**

**Wheels:** 2  
**Power Source:** 100-240V adapter included, 6 D batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Audible / Mute Sound Option, Everlock®

Model	Description
2500KG	Scale in kilograms only
2500KL-BT	2500KL with built-in Pelstar® wireless technology
2500KG-BT	2500KG with built-in Pelstar® wireless technology
SS-2500KL	2 year extended warranty
CHAIRACC	Optional fold away seat
ADPT30	Adapter included, item # for replacement or spare
B2500RAMP	Secondary ramp for 2500KL and 2500KG
BTPRINTERKIT-01	Wireless thermal printer kit



Facilities looking for a scale that can be used with wheelchairs and takes advantage of a built-in chair, should consider the 2500CKL, which features both a ramp and a folding patient seat.

**Platform Capacity:** 1000 lb / 454 kg  
**Seat Capacity:** 419 lb / 190 kg  
**Resolution:** 0.2 lb / 0.1 kg  
**Platform Size (w x d x h):** 29 ⅜" x 28 ¼" x 2 ⅝"  
746 mm x 718 mm x 60 mm  
**Seat Size (w x d x h):** 18" x 15" x 18 ½"  
457 mm x 381 mm x 470 mm  
**Display:** Digit Size: ⅜" / 10 mm, Screen Size: 2 ¾" x 1 ½" High-Contrast Color TFT-LCD Screen  
**Product Footprint (w x d x h):** 40 ¼" x 34 ½" x 45 ½"  
1022 mm x 876 mm x 1156 mm  
**Options & Accessories:**


**Connectivity:** USB, optional Pelstar® wireless technology  
**Wheels:** 2  
**Power Source:** 100-240V adapter included, 6 D batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Audible / Mute Sound Option, Everlock®

Model	Description
2500CKG	Scale in kilograms only
SS-2500CKL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
B2500RAMP	Secondary ramp for 2500CKL
BTPRINTERKIT-01	Wireless thermal printer kit



2650KL

DIGITAL WHEELCHAIR DUAL RAMP SCALE

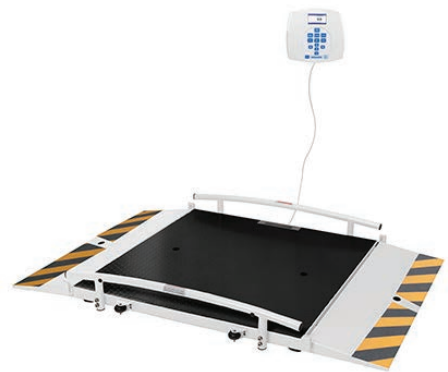

Facilities frequently weighing patients in wheelchairs will appreciate the 2650KL, which offers a remote display head and a medium-sized platform. The dual ramps make it easy to wheel patients on and off the platform. Like all wheelchair scales, the 2650KL can also be used for walk-on patients.

- Capacity:** 1000 lb / 454 kg

**Resolution:** 0.2 lb / 0.1 kg

**Platform Size (w x d x h):**  
31 ½" x 31 ½" x 2"  
800 mm x 800 mm x 51 mm

**Display:** Digit Size: ¾" / 10 mm,  
Screen Size: 2 ¾" x 1 ½" High-Contrast  
Color TFT-LCD Screen

**Product Footprint (w x d x h):**  
51 ¾" x 33 ⅝" x 5 ¼"  
1314 mm x 854 mm x 133 mm

**Options & Accessories:**
- Connectivity:** USB, optional Pelstar® wireless technology

**Wheels:** 2

**Power Source:** 100-240V adapter included, 6 D batteries (not included)


**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Everlock®, Audible / Mute Sound Option

Model	Description
2650KG	Scale in kilograms only
2650KL-BT	2650KL with built-in Pelstar® wireless technology
2650KG-BT	2650KG with built-in Pelstar® wireless technology
SS-2650KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit



2700KL

DIGITAL WHEELCHAIR DUAL RAMP SCALE WITH EXTRA LARGE PLATFORM

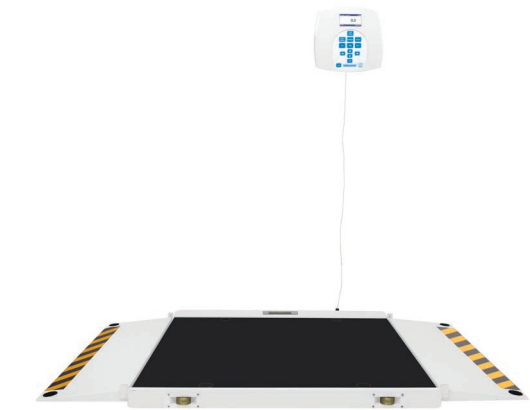

Featuring the largest platform in the wheelchair scale collection, the 2700KL can easily accommodate even the widest wheelchair. It is a good choice in facilities where the scale will remain in one location.

- Capacity:** 1000 lb / 454 kg

**Resolution:** 0.2 lb / 0.1 kg

**Platform Size (w x d x h):**  
43" x 42" x 2 ½"  
1092 mm x 1067 mm x 64 mm

**Display:** Digit Size: ¾" / 10 mm,  
Screen Size: 2 ¾" x 1 ½" High-Contrast  
Color TFT-LCD Screen

**Product Footprint (w x d x h):**  
69 ⅞" x 51 ½" x 3 ⅜"  
1756 mm x 1308 mm x 86 mm

**Options & Accessories:**
- Connectivity:** USB, optional Pelstar® wireless technology

**Wheels:** 2

**Power Source:** 100-240V adapter included, 6 D batteries (not included)

**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Everlock®, Audible / Mute Sound Option

Model	Description
2700KG	Scale in kilograms only
2700KL-BT	2700KL with built-in Pelstar® wireless technology
2700KG-BT	2700KG with built-in Pelstar® wireless technology
SS-2700KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit



2900KL-AM

ANTIMICROBIAL WALL-MOUNTED FOLD-UP WHEELCHAIR SCALE



- Power Source:** 100-240V adapter included, 6 D batteries (not included)
- Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Everlock®, Audible / Mute Sound Option
- Connectivity:** USB, optional Pelstar® wireless technology

2920KL-AM

ANTIMICROBIAL WALL-MOUNTED FOLD-UP WHEELCHAIR SCALE, ADA COMPLIANT



Learn More



32" x 40" low-profile platform, 2" edge protection, and 1:8 running slope complies with MDE standards

The 2900KL-AM is the perfect scale for facilities that need to save space. The low profile scale extends off the wall 3.5" when folded up and has a large easy-to-lift platform that accommodates most wheelchairs.

**EMR Capable:** Scale interfaces with a computer, EMR software or other electronic device via Health o meter® Professional Connectivity Solutions.

**Easy to Lift Large Platform:** Large easy-to-lift platform accommodates most wheelchairs and provides a comfortable weighing experience.

**Safety Features:** Anti-slip caution tape, antimicrobial protection, raised platform edge, recessed handle and no pinch points provide both the clinician and patient confidence during the weighing experience.

Capacity	Resolution	Platform Size (w x d x h)
1000 lb / 454 kg	0.2 lb / 0.1 kg	29 ½" x 32 ½" x 2 ¾" 749 mm x 819 mm x 70 mm
Display	Product Footprint (w x d x h)	
Digit Size: ¾" / 10 mm, Screen Size: 2 ¾" x 1 ½" High-Contrast Color TFT-LCD Screen	58" x 37" x 37" (unfolded) 1473 mm x 940 mm x 940	

Options & Accessories:

Model	Description
2900KG-AM	Scale in kilograms only
2900KL-AM-BT	2900KL- AM with built-in Pelstar® wireless technology
2900KG-AM-BT	2900KG-AM with built-in Pelstar® wireless technology
2920KL-AM	ADA compliant scale
2920KG-AM	ADA compliant scale in kilograms only
SS-2900KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit



The 2920KL-AM is an ideal solution for facilities that need to meet accessibility standards while conserving space. The low profile scale extends off the wall 4.7" when folded up and has a large easy-to-lift platform that accommodates most wheelchairs.

- Capacity:** 1000 lb / 454 kg

**Resolution:** 0.2 lb / 0.1 kg

**Platform Size (w x d x h):**  
40" x 32" x 1 ⅝"  
1016 mm x 819 mm x 38 mm

**Display:** Digit Size: ¾" / 10 mm,  
Screen Size: 2 ¾" x 1 ½" High-Contrast  
Color TFT-LCD Screen

**Product Footprint (w x d x h):**  
77 ½" x 36 ½" x 37" (unfolded)  
1968 mm x 927 mm x 940 mm

**Options & Accessories:**
- Connectivity:** USB, optional Pelstar® wireless technology

**Power Source:** 100-240V adapter included, 6 D batteries (not included)

**Functions:** LB / KG Conversion, LB / KG Lock Out, BMI, Zero, Tare, Pre-tare, Hold / Release, Recall, Reweigh, Auto Zero, Auto Off, Time / Date, Variable Auto Off Time, Everlock®, Audible / Mute Sound Option

Model	Description
2920KG-AM	Scale in kilograms only
SS-2920KL-AM	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit





# Chair Scales

## 2595KL

DIGITAL CHAIR SCALE



This best-selling chair scale has a comfortable seat so patients can sit safely while being weighed. The moveable armrests and footrests allow for easy access on and off the scale. A great scale for use in Long-Term Care centers, rehab facilities, and senior living homes.

**Capacity:** 600 lb / 272 kg  
**Resolution:** 0.2 lb / 0.1 kg  
**Seat Size** (w x d x h):  
18 ¼" (w) x 15 ⅝" (d) x 16 ½"  
464 mm x 397 mm x 419 mm  
**Display:** 1 ⅞" / 29 mm LCD  
**Connectivity:** optional Pelstar®  
wireless technology  
**Product Footprint** (w x d x h):  
23 ½" x 37 ½" x 39 ⅞"  
596 mm x 952 mm x 1012 mm

**Options and Accessories:**

Model	Description
2595KG	Scale in kilograms only
SS-2595KL	2 year extended warranty
ADPT30	Adapter included, item # for replacement or spare



**Wheels:** 2  
**Power Source:** 100-240V adapter included, 6 C-cell batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, Zero, Pre-tare, Hold / Release, Recall, Reweigh, Tare, Auto Off, Variable Auto Off Time, Enabling Unit View with Unit of Measure Locked, Everlock®

## 594KL

DIGITAL CHAIR SCALE



The 594KL has moveable arm and foot rests to allow for easy patient access on and off this economical chair scale. Four casters make portability easy for bringing the scale to the patient's location.

**Capacity:** 600 lb / 270 kg  
**Resolution:** 0.2 lb / 0.1 kg  
**Seat Size** (w x d):  
18" x 18 c"  
457 mm x 470 mm  
**Display:** 1" / 25 mm LCD  
**Product Footprint** (w x d x h):  
23 ¼" x 40 ⅞" x 37 ⅞"  
590 mm x 1020 mm x 960 mm

**Options and Accessories:**

Model	Description
594KL-C	Scale, Canada only
SS-594KL	2 year extended warranty
69-00001	Adapter included, item # for replacement or spare

**Wheels:** 4 casters  
**Power Source:** 100-240V adapter included, 6 AA batteries (not included)  
**Functions:** LB / KG Conversion, LB / KG Lock Out, BSA, BMI, Zero, Tare, Hold / Release, Auto Zero, Auto Off

# Pediatric Scales and Carts

★TOP SELLER

## 553KL

DIGITAL PEDIATRIC TRAY SCALE  
WITH EXTRA-WIDE TRAY



Optional Digital and Mechanical Height Rods Available



The 553KL is Health o meter® Professional's most popular infant scale and perfect for newborns, infants and toddlers. The scale operates on batteries or AC power so you can use it anywhere. Practices benefit from one-button operation that makes it easy to use while ensuring that the baby remains safely on the scale.

**Measuring Tape Included:** Easy-to-read measuring tape allows user to measure both weight and height on the same device.

**Motion Detection:** Scale weighs accurately despite movement of the child on the scale, and locks onto the measurement.

Capacity	Resolution	Tray Size (w x d x h)
----------	------------	-----------------------

44 lb / 20 kg	0.5 oz / 10 g	24 ⅞" x 13" x 1 ⅝" 632 mm x 330 mm x 35 mm
---------------	---------------	---

Measuring Tape	Display	Product Footprint (w x d x h):
----------------	---------	--------------------------------

0" - 24" 0 cm - 58 cm	1" / 25 mm LCD	24 ⅞" x 16 ½" x 3 ⅞" 632 mm x 419 mm x 98 mm
--------------------------	----------------	---

**Options & Accessories:**

Model	Description
553KLCT	Scale with rolling cart included
553KG	Scale in kilograms only
553KL-HR	553KL with attachable mechanical height rod
553KG-HR	553KG with attachable mechanical height rod
553KL-EHR	553KL with attachable digital height rod
553KG-EHR	553KG with attachable digital height rod
SS-553KL	2 year extended warranty
553HR	Attachable mechanical height rod for 553 series scales
553EHR	Attachable digital height rod for 553 series scales
2210CART	Rolling cart for 553KL & 553KG
553CASE	Carrying case for 553KL & 553KG
ADPT31	Adapter for 553KL & 553KG
BTTPRINTERKIT-01	Wireless thermal printer kit



**Power Source:** 6 AA batteries included, 100-240V adapter (not included)  
**Functions:** LB / KG Conversion, LB / KG Lockout, Zero, Tare, Pre-Tare, Hold / Release, Reweigh, Recall, Auto Zero, Auto Off, Motion Detection, Variable Auto-Off Time, Enabling Unit View with Unit of Measure Locked, Auto-Weight Hold for 30 to 60 Seconds, Everlock®  
**Connectivity:** optional Pelstar® wireless technology



2210CART

2210KL-AM

ANTIMICROBIAL HIGH RESOLUTION  
DIGITAL NEONATAL / PEDIATRIC TRAY  
SCALE



2210KL-AM  
with 2210EHR



2210CART-AC

2210KL4-AM

ANTIMICROBIAL HIGH RESOLUTION  
DIGITAL NEONATAL/PEDIATRIC SCALE  
WITH 4-SIDED TRAY



2210KL4-AM

The 2210KL-AM neonatal scale was designed for use in the NICU, delivery room, nursery and pediatric practices where a high degree of accuracy is critical to ensure that babies are gaining weight. A mobile cart designed for the 2210KL-AM series makes it easy to move the scale between locations.

**Capacity:** 45 lb / 20 kg / 20,000 gm  
**Resolution:** 1 g (to 10,000 g);  
5 g (10,001 to 20,000 g)

**Tray Size** (w x d x h):  
24" x 15 ¼" x 2 ¾"  
610 mm x 388 mm x 71 mm

**Display:** 5 ¼" x 1 ⅛" (133 mm x 35 mm),  
4 levels of brightness

**Connectivity:** via USB, optional Pelstar®  
wireless technology

**Product Footprint** (w x d x h):  
24 ¾" x 17 ¾" x 6"  
630 mm x 450 mm x 153 mm

Options & Accessories

Model	Description
2210KG-AM	Scale in kilograms only
2210KL-AM-BT	2210KL-AM with built-in Pelstar® wireless technology
2210KG-AM-BT	2210KG-AM with built-in Pelstar® wireless technology
2210KL-AM-C	Scale, Canada only
2210KG-AM-C	Scale in kilograms only, Canada only
2210KL-AM-BT-C	2210KL-AM with built-in Pelstar® wireless technology, Canada only
2210KG-AM-BT-C	2210KG-AM with built-in Pelstar® wireless technology, Canada only
2210EHR	Attachable digital height rod for 2210KL-AM (2-sided tray scale)
SS-2210KL	2 year extended warranty
2210CART-AC	Acute care cart for 2210KL-AM series of scales
ADPT30	Adapter included, item # for replacement or spare
BTTPRINTERKIT-01	Wireless thermal printer kit

The 2210KL4-AM neonatal scale comes with a 4-sided weighing tray, which provides an additional level of safety and protection for the baby during the weighing process.

**Capacity:** 45 lb / 20 kg / 20,000 gm  
**Resolution:** 1 g (to 10,000 g);  
5 g (10,001 to 20,000 g)

**Tray Size** (w x d x h):  
26 ½" x 15 ¼" x 2 ¾"  
674 mm x 388 mm x 71 mm

**Display:** 5 ¼" x 1 ⅛" (133 mm x 35 mm),  
4 levels of brightness

**Connectivity:** via USB, optional Pelstar®  
wireless technology

**Product Footprint** (w x d x h):  
26 ½" x 17 ¾" x 6"  
674 mm x 450 mm x 153 mm

Options & Accessories:

Model	Description
2210KG4-AM	Scale in kilograms only
2210KL4-AM-BT	2210KL4-AM with built-in Pelstar® wireless technology
2210KG4-AM-BT	2210KG4-AM with built-in Pelstar® wireless technology
2210KL4-AM-C	Scale, Canada only
2210KG4-AM-C	Scale in kilograms only, Canada only
2210KL4-AM-BT-C	2210KL4-AM with built-in Pelstar® wireless technology, Canada only
2210KG4-AM-BT-C	2210KG4-AM with built-in Pelstar® wireless technology, Canada only
SS-2210KL	2 year extended warranty
2210CART-AC	Acute care cart for 2210KL-AM series of scales
ADPT30	Adapter included, item # for replacement or spare
BTTPRINTERKIT-01	Wireless thermal printer kit

**Measuring Tape:** 0" - 24" / 0 cm - 60 cm  
**Power Source:** 6 D batteries included,  
100-240V adapter included

**Functions:** LB / KG / GM Conversion,  
LB / KG / GM Lock Out, Zero, Tare,  
Pre-tare, Reweigh, Recall, Auto Zero,  
Auto Off, Motion Detection, Everlock®,  
Hold / Release, Variable Auto Off Time,  
Audible / Mute Sound Option, Disable  
/ Enable Backlight for Display, Auto  
Weight Hold for 30 to 60 seconds, Breast  
Milk Intake

522KL

DIGITAL PEDIATRIC TRAY SCALE



This easy-to-use scale is perfect for busy pediatric offices or clinics where infants are seen regularly. The tray has a convenient measuring tape to determine the baby's length during the weighing session.

**Capacity:** 50 lb / 23 kg  
**Resolution:** 0 - 20 lb / 0.2 oz; 20 - 50 lb /  
0.5 oz; 0 - 9 kg / 5 g; 9 - 23 kg / 10 g

**Tray Size** (w x d x h):  
24 ⅞" x 14 ½" x 2 ⅝"  
613 mm x 368 mm x 67 mm

**Measuring Tape:** 0" - 23" / 0 cm - 58 cm

**Connectivity:** optional Pelstar®  
wireless technology

**Product Footprint** (w x d x h):  
24 ⅞" x 21" x 23 ⅞"

613 mm x 533 mm x 587 mm

Options & Accessories:

Model	Description
522KG	Scale in kilograms only
522KL-HR	522KL with attachable mechanical height rod
522KG-HR	522KG with attachable mechanical height rod
522KL-EHR	522KL with attachable digital height rod
522KG-EHR	522KG with attachable digital height rod
SS-522KL	2 year extended warranty
522HR	Attachable mechanical height rod for 522 series scales
522EHR	Attachable digital height rod for 522 series scales
ADPT31	Adapter for 522 series scales
BTTPRINTERKIT-01	Wireless thermal printer kit



524KL

DIGITAL PEDIATRIC SEAT SCALE



The deep infant seat on the 524KL reduces the movement of young children to speed up the weighing process. Weight is captured easily while the seat supports the child's head and body.

**Capacity:** 50 lb / 23 kg  
**Resolution:** 0 - 20 lb / 0.2 oz; 20 - 50 lb /  
0.5 oz; 0 - 9 kg / 5 g; 9 - 23 kg / 10 g

**Seat Size** (w x d x h):  
24 ¼" x 14 ¼" x 18 ½"  
616 mm x 362 mm x 470 mm

**Connectivity:** optional Pelstar® wireless  
technology

**Product Footprint** (w x d x h):  
24 ½" x 21" x 23 ⅞"

622 mm x 533 mm x 587 mm

Options & Accessories:

Model	Description
524KG	Scale in kilograms only
SS-524KL	2 year extended warranty
ADPT31	Adapter for 524KL & 524KG





565KL

HIGH RESOLUTION DIGITAL  
BABY SCALE



565KL  
with 565EHR

The 565KL was designed for use in critical settings such as the NICU, nursery, delivery room, and pediatric practices. Its exceptional accuracy is especially crucial when caring for the tiniest patients.

- Capacity:** 44 lb / 20 kg

**Resolution:** 0 – 44 lb x 0.1 oz / 0 – 6,000 g x 1g / 6,001 g – 20,000 g x 2 g

**Tray Size (w x d x h):**  
24 7⁄8" x 9 7⁄8" x 2 3⁄4"  
632 mm x 250 mm x 70 mm

**Connectivity:** optional Pelstar® wireless technology

**Product Footprint (w x d x h):**  
26 1⁄4" x 13" x 5 1⁄4"  
(667 mm x 330 mm x 134 mm)

**Options & Accessories:**
- Display:** Digit Size: 3⁄4" / 19 mm, Screen Size: 3" x 1" LCD display

**Measuring Tape:** 0" - 24 5⁄8" / 0 cm - 62.5 cm

**Power Source:**  
100-240V Adapter included, 6 AA batteries not included

**Functions:** LB / KG Conversion, LB / KG Lock Out, Zero, Tare, Pre-tare, Reweigh, Recall, Auto Zero, Auto Off, Motion Detection, Everlock®, Hold / Release, Variable Auto Off Time, Audible / Mute Sound Option, Auto Weight Hold for up to 5 Minutes

Model	Description
565KG	Scale in kilograms only
565KL-C	Scale, Canada only
565KG-C	Scale in kilograms only, Canada only
565EHR	Attachable digital height rod for 565 series scales
2210CART	Rolling cart for 565KL & 565KG
SS-565KL	2 year extended warranty
69-00001	Adapter included, item # for replacement or spare
BTPRINTERKIT-01	Wireless thermal printer kit

Pediatric Carts

2210CART-AC

ACUTE CARE CART FOR 2210KL-AM / 2210KL4-AM  
SERIES OF SCALES



- Features:** 2 easy grip handles and customized plates for secure scale placement of the 2210KL-AM series of scales
- Storage Shelf (w x d x h):**  
21 1⁄4" x 13" x 3"  
540 mm x 330 mm x 76 mm
- Capacity:** 77 lb / 35 kg
- Cart Dimensions (w x d x h):**  
30" x 22" x 38"  
762 mm x 59 mm x 965 mm
- Wheels:** 4 large anti-static casters with locks

2210CART

ROLLING SCALE / MEDICAL DEVICE CART FOR USE WITH  
549KL, 550KL, 553KL, 565KL, 2210KL-AM, 2210KL4-AM



- Features:** Lockable cabinet and storage shelf.
- Cabinet Dimensions (w x d x h):**  
23" x 16 3⁄4" x 17"  
584 mm x 427 mm x 432 mm
- Capacity:** 200 lb / 90 kg
- Cart Dimensions (w x d x h):**  
26 1⁄2" x 18" x 34 1⁄8"  
673 mm x 457 mm x 867 mm
- Wheels:** 4 casters (2 with locks)



Pediatric Scales Designed for  
Safety & Workflow

Specialty Scales

★TOP SELLER

PTS-1000KG

Patient Transfer Scale, Kilograms Only

The PTS-1000KG a transfer board with a built-in medical scale, is designed to significantly reduce treatment times in stroke protocols or other time-sensitive medical conditions where precise weight measurement of immobile patients is crucial. Lightweight and portable, the PTS integrates seamlessly into trauma bay and other clinical workflows, providing a rapid, accurate weight at the point of care, without delaying treatment or requiring the patient to be transferred off the stretcher. Health o meter® Professional Scales warrants that the Patient Transfer Scale will be accurate to within 500 grams so long as the scale is properly used and maintained.

**Designed by a Clinician for Clinicians:** Weigh critical or immobile patients instantly and accurately helping shorten door-to-needle times.



**Power Source:** 3.7V rechargeable battery  
**Functions:** Zero, Hold / Release, Auto Off

Always use a Slide Sheet when  
using the Patient Transfer Scale



Check out the Single Patient Use  
and Reusable SlideMate™

Capacity	Resolution	Weighing Mat Size (w x d):
250 kg	0.5 kg	71" x 21" 1805 mm x 533 mm
Measuring Tape	Display	Product Footprint (w x d x h):
100 cm -180 cm height guide on mat	3" x 1" Backlight LCD Screen	71" x 27 1⁄2" x 1" 1805 mm x 700 mm x 30 mm

Options & Accessories:

Model	Description
PTS-1000KL	Scale, lbs and kilograms
PTS-1000KL-C	Scale, lbs and kilograms, Canada only
PTS-1000KG-C	KG only scale, Canada only
SS-PTS-1000KL	2 year extended warranty
69-00012	Battery charger included, item # for replacement or spare



BCS-G7

BODY COMPOSITION SCALE



**Power Source:** 110-220V power adapter (included)  
**Functions:** 30 second testing start to finish  
**Full Report Provides:** Weight, Body fat %, Body water %, Muscle mass %, Fat Free Mass, Ideal Weight Range, Segmental Fat Distribution Analysis, Total Body Fat, Total Body Water, BMI – Body Mass Index, BMR – Basal Metabolic Rate, REE – Resting Energy Expenditure, Active Metabolic Rate, Visceral Fat

This body composition scale features an industry-leading 1000 lb capacity, that provides the flexibility and accuracy to weigh any patient. Few pretest requirements that the analysis can be performed at any time of day and under almost any circumstances, resulting in the greatest flexibility of any body composition analyzer.

**15 Second Measurement Time:** On-screen display delivers rapid access to weight, body fat %, body water %, lean mass, BMR & BMI. A full illustrated printout can be generated in under a minute.

**Large Platform:** The spacious platform allows patients to find a comfortable and secure stance.

**Few Pretest Restrictions:** Analysis can be performed at any time of day and under almost any circumstances, resulting in the greatest flexibility of any body composition analyzer.

**High Capacity:** The scale features an industry-leading 1000 lb capacity that provides the flexibility and accuracy to weigh any patient.

Capacity	Resolution	Platform Size (w x d x h)
1000 lb / 453.5 kg	0.1 lb / 0.1 kg	14" x 16" x 2 7⁄8" 356 mm x 406 mm x 48 mm

Display	Product Footprint (w x d x h):
6-window display (.5" digits) includes Weight, % Body Fat, % Body Water, % Muscle Mass, BMI, REE	20" x 18" x 41" 508 mm x 457 mm x 1041 mm



3400KL / 3401KL

DIGITAL WET DIAPER / LAP SPONGE / ORGAN SCALE



**Power Source:** 6V rechargeable battery included, 100 - 240V AC battery charger included  
**Functions:** Zero, Tare, Auto Zero, Auto Off, g / oz Conversion, Low Battery Indicator, Backlit Display

The solid, industrial design of the 3400KL includes a stable, stainless steel base and removable, stainless steel weighing pan (3401KL comes with a flat tray) that can be autoclaved for sterilization. Powered by a long-life rechargeable battery and backlit for any lighting level, this portable scale can be used in any environment, from OR to ED to NICU and beyond.

<b>Capacity:</b> 176 oz / 5,000 g	<b>Tray Size (3401KL)</b> (w x d x h): 7 ¾" x 9 ¾" x 5⁄8"
<b>Resolution:</b> 0.05 oz / 1 g	197 mm x 248 mm x 16 mm
<b>Base Size</b> (w x d x h): 11 ½" x 9 7⁄8" x 5 ¼"	<b>Display:</b> 5⁄8" / 16 mm LCD
292 mm x 250 mm x 133 mm	<b>Product Footprint</b> (w x d x h): 10 7⁄8" x 13" x 8" (3400KL)
<b>Pan Size (3400KL)</b> (w x d x h): 10" x 9 7⁄8" x 2 ½"	276 mm x 330 mm x 203 mm
254 mm x 250 mm x 62 mm	10" x 11 ¾" x 5 3⁄8" (3401KL)
	254 mm x 298 mm x 137 mm

Options & Accessories:

Model	Description
3400KG	Scale with pan, kilograms only
3401KL	Scale with tray
3400KG-C	Scale with pan, kilograms only, Canada only
SS-3400KL	2 year extended warranty
SS-3400KG	2 year extended warranty
SS-3400KG-C	2 year extended warranty
SS-3401KL	2 year extended warranty
3400PAN	3 pack of stainless steel pans
3400TRAY	3 pack of stainless steel trays
3400ADPT	Battery charger for 3400KL, 3400KG, & 3401KL
20-00018	Battery charger for 3400KG-C



Height Measuring

201HR

WALL-MOUNTED TELESCOPIC METAL HEIGHT ROD



★ TOP SELLER

Sleek and robust measuring device requires minimal wall space and has a foldable headpiece. Designed for long life and to withstand dozens of uses every day.

**Measurement Range:** 2" - 90" / 6 cm - 230 cm  
**Graduation:** 1/16" / 1 mm  
**Product Footprint** (w x d x h): 2 ¼" x ¾" x 54", 57 mm x 19 mm x 1372 mm

205HR

HIGH-STRENGTH WALL-MOUNTED HEIGHT ROD



Manufactured with a durable high-strength plastic polymer, this sturdy stadiometer resists bending and ensures that the headpiece slides smoothly. Featuring a double sided measuring tape, it offers a choice of centimeters only or feet and inches.

**Measurement Range:** 55.5" / 141 cm  
**Measuring Tape:** 0" - 90" / 0 cm - 230 cm  
**Graduation:** 1/16" / 1 mm  
**Product Footprint** (w x d x h): 2 5⁄8" x ¾" x 55", 66 mm x 19 mm x 1397 mm



280HR

INFANT MEASURING MAT



The space-saving mat is soft, lightweight, foldable, and easy to store. The fixed headpiece and sliding foot positioner allow for precise measurement of infants or toddlers while lying down.

**Measurement Range:** 3.9" - 43.3" / 10 cm - 110 cm  
**Graduation:** 1/4" / 5 mm  
**Product Footprint** (l x w): 55" x 12", 1397 mm x 305 mm

EBABYHR

DIGITAL PEDIATRIC MEASURING ROD



The eBABYHR makes measuring infant length fast and easy with clear digital measurements, large head and foot paddles, and convenient Hold function.

**Measurement Range:** 13 ¾" - 31 ½" / 35 cm - 80 cm  
**Resolution:** 1/16" / 0.1 cm  
**Product Footprint** (w x d x h): 24 ¾" x 11 ½" x 2 ¾", 629 mm x 292 mm x 69 mm



# Scale Service and Maintenance

## Battery Life & Power Adapters

- Scales have the option to be powered by batteries or an AC adapter. Most scales also feature an Auto Off function that automatically powers off the scale after a period of time to preserve battery life.

## Zeroing

- Some organizations direct staff to zero a scale before each patient weighing. Health o meter® Professional scales automatically zero when the scale is powered on but most also feature a Zero button to manually zero the scale.

## Replacement Parts

- Health o meter® Professional scales are designed to withstand frequent use in any medical setting. However, should a part of the scale become damaged replacement parts are available.

## Revit Files

- Technical files for the most popular scales are available for download:  
[www.hom scales.com/product-support/revit-files/](http://www.hom scales.com/product-support/revit-files/)

## Calibration

- Health o meter® Professional scales are factory calibrated. We do understand that healthcare facilities have their own calibration standards and when necessary all of our scales can be easily calibrated. Calibration procedures can be found at [www.hom scales.com/technical-support/calibration/](http://www.hom scales.com/technical-support/calibration/)

## Calibot


- For large organizations Health o meter® Professional offers an innovative, exclusive device to make calibration weigheasier. The Automated Calibrator is a robust innovative calibration system that allows biomedical engineers and other medical technicians to calibrate medical patient scales quickly and safely without having to transport and stack physical calibration weights on the scale. Eliminating the cumbersome physical calibration weights from the calibration process reduces the time needed for routine scale maintenance and eliminates the safety risks associated with handling heavy dead weights.



# Customer Support

We have over 170 skilled professionals based in the United States and Canada who manufacture, sell, and support our products. Our global headquarters in McCook, IL is home to our highly trained Customer Service and Technical Support teams along with our engineering, logistics and management team. On-site sales and support are provided in partnership with MedTech/MedCare, an industry-recognized leader in medical sales throughout the US and Canada.

- Direct 800# for Technical Support – Representatives available to answer any functionality questions. All calls answered in less than 2 rings.
- In-House Service Specialists – With 100% focus on medical scales, each team member is a product specialist and can assist quickly and easily.
- On-Site Support – Local representatives can assist with product guidance and recommendations, staff training, and facility equipment evaluations.

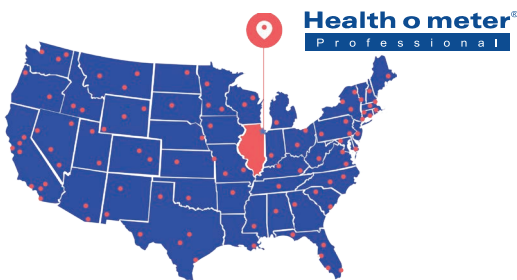

## Where to Purchase

Health o meter® Professional scales are available through all major medical supply distributors as well as multiple Group Purchasing Organizations. Sales guidance is also provided through local representatives of MedTech/MedCare. To find your local rep visit [www.mtmhealth.com/find-a-local-representative/](http://www.mtmhealth.com/find-a-local-representative/)

## MTMC

### About MedTech/MedCare

MedTech/MedCare is our trusted independent sales team, bringing decades of experience and a deep understanding of the healthcare industry. With a strong commitment to service and innovation, they partner closely with healthcare providers, systems, and distributors to deliver tailored solutions that support better outcomes and operational efficiency. Learn more about their legacy and mission at [www.mtmhealth.com](http://www.mtmhealth.com)

## Warranty

Pelstar, LLC/Health o meter® Professional scales are warranted from date of purchase against defects of materials or in workmanship for a period of time (see product instruction manual for the warranty period). If product fails to function properly, return the product, freight prepaid and properly packed to Pelstar. If manufacturer determines that a defect of material or workmanship exists, a customer's sole remedy will be replacement of scale at no charge. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. All parts including replaced parts are covered only for the original warranty period.

**Ask Us About AIA Accredited Education for your Firm!**



# Health o meter®

P r o f e s s i o n a l

## Contact Information

Main Number: 1 (708) 377-0600

Customer Service: 1 (800) 815-6615

Technical Support: 1 (800) 638-3722

Fax: 1 (708) 377-0601 or 1 (800) 654-7330

Email: [homprocs@homscales.com](mailto:homprocs@homscales.com)

Corporate Headquarters:

Pelstar LLC

9500 West 55th St.

McCook, IL 60525-7110 USA

[www.homscales.com](http://www.homscales.com)

Health o meter® Professional products are manufactured, designed and owned by Pelstar LLC. Health o meter® is a registered trademark of Sunbeam, Inc. used under license. ©2025 Pelstar® LLC – SCALE\_PLANNING\_GUIDE\_20251203