

## THE 2622KL

# DIGITAL WHEELCHAIR AND STAND-ON SCALE WITH DUAL RAMPS AND EXTRA-LARGE PLATFORM, CERTIFIED FOR ADA COMPLIANCE

The all-digital 2622KL is engineered with enhanced safety features to ensure a comfortable weighing experience for patients standing on the platform or for those who use wheelchairs. Certified for Americans with Disabilities Act (ADA) compliance, the dual ramps adhere to a 1:8 slope ratio and the handrail height is within the required 34" -38" range. The ramps and extra-large 32" x 40" platform feature a 2" edge protection to prevent roll-offs, and anti-slip caution tape on the ramps provide enhanced safety. The platform folds for storage, and built-in casters are included for added convenience and mobility. The scale has a weight capacity of up to 1,000 pounds (454 kilograms). As an all-in-one scale, the 2622KL meets the U.S. Access Board Standards for both stand-on and wheelchair scales, and is ideal for doctor's offices, hospitals, clinics, and acute or long-term care facilities. Additional features include an antimicrobial keypad, a wall-mounted display head, and a large non-skid platform.



The scale features a high-contrast color screen with user-friendly functions like BMI, hold/release, reweigh, recall, zero, KG/LB, tare, and pre-tare. The Everlock® function permanently locks the measuring unit (LB or KG). With a USB interface, it seamlessly integrates with EMR systems and connects to devices via Health o meter® Professional Connectivity Solutions.

The 2622KL is designed to meet the US Access Board's Standards. Medical facilities operated by state and local governments, also known as public entities, must comply with the MDE Standards for Accessible Medical Diagnostic Equipment, which are enforced by the Department of Justice.

- 32" wide x 40" deep platform
- 2" edge protection on platform
- 2" edge protection on ramps
- 1:8 running slope
- Slip-resistant platform
- Handrail support along the length of the platform
- 34" - 38" max. handrail height



**Health o meter®**  
Professional

Health o meter® Professional 2622KL

Reviewed for  
**ADA compliance and  
Platinum certified by Corada**

## THERE'S STRENGTH IN OUR NUMBERS.

When you add up the numbers on our scale, it's not surprising that the 2622KL provides a safe, accessible weighing experience.



**100%**  
CERTIFIED FOR  
ADA COMPLIANCE

**2"** EDGE PROTECTION

**32" X 40"** LOW  
PROFILE PLATFORM

**37.5"** HANDRAIL  
HEIGHT

**1:8** RUNNING  
SLOPE

**#1** TRUSTED  
BRAND

weigh**easier**®

# TECHNICAL SPECIFICATIONS

## Capacity

1000 lb / 454 kg

## Resolution

0.2 lb / 0.1 kg

## EMR connectivity

Via USB, optional Pelstar® wireless technology

## Display

Digit Size: 3/8" / 10 mm  
Screen Size: 2 3/4" x 1 1/2" High-Contrast  
Color TFT-LCD Screen

## Power Source

100-240V adapter (ADPT30) included, 6 D batteries not included

## Functions

LB/KG conversion  
LB/KG lock out  
BMI  
Zero  
Tare  
Pre-tare  
Hold/release  
Reweigh  
Recall  
Auto zero  
Auto off  
Time/date  
Variable auto off time  
Audible/mute sound option  
Parent/child weighing function  
Everlock®

## Platform Size (w x d x h)

40" x 32 1/4" x 1 7/8"  
(1016 mm x 819 mm x 48 mm)

## Ramp Size

32 1/4" x 16"  
(819 mm x 406 mm)

## Edge Guard Size (w x d x h)

40" x 1/8" x 2"  
(1016 mm x 3 mm x 51 mm)

## Product Footprint (w x d x h)

(In Operating Position)  
71 1/2" x 45 1/8" x 39 1/2"  
(1816 mm x 1147 mm x 1003 mm)

## Product Footprint (w x d x h)

(Folded)  
71 1/2" x 9 1/2" x 44"  
(1816 mm x 241 mm x 1118 mm)

## Product Weight

126 lb / 57 kg

## Handrail Height

37.5" / 952 mm

## Portability

4 wheels (2 casters)  
for ease of movement

## Warranty


2-year limited manufacturer warranty

# ORDERING INFORMATION

MODEL	DESCRIPTION
2622KL	Digital Wheelchair and Stand-On Scale with Dual Ramps and Extra-Large Platform, ADA Compliant, LB & KG*
2622KG	Digital Wheelchair and Stand-On Scale with Dual Ramps and Extra-Large Platform, ADA Compliant, KG only*
43-0070	Optional Transport Kit for 2622 Series Display Head
SS-2622KL	ScaleSurance® 2 year extended warranty for 2622 series*
ADPT30	Adapter included, item # for replacement or spare
PLUGSET30	Plug adapters for use outside North America
BTPRINTERKIT-01	Wireless thermal printer kit

Scale also includes connectivity capabilities. To learn about connectivity options, visit [hom scales.com/ConnectivitySolutions](http://hom scales.com/ConnectivitySolutions)

To learn more, call **1-800-815-6615**  
or visit **hom scales.com/2622KL**

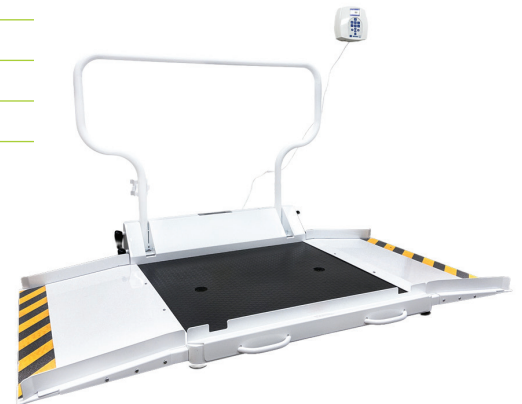

**Health o meter®**  
Professional



\*Not available in all countries.